WATERS ANALYTICAL STANDARDS AND REAGENTS

Precision Defined.



THE SCIENCE OF WHAT'S POSSIBLE.™

Waters commitment

Waters understands that the quality of the standards and reagents you use directly correlates to the quality of your results and final products. That is why we have made analytical standards and reagents our business. Our seal symbolizes Waters commitment to you.

Precise Formulation

In the field of analytical measurement, the ability to repeatedly prepare reference materials identically over long periods of time is critical for data comparability and defensibility.

Absolute Traceability

In order to be true, the properties of the measurement must be directly linked to a source reference material through an unambiguous, unbroken and fully documented chain of comparisons.

Only from Waters

A keen understanding of reagents and certified reference materials, together with the ability to prepare, test and deliver them as certified, ready-to-use products is unique to one analytical instrumentation manufacturer, Waters.



Waters Standards and Reagents Products

Standards for System Performance
MS Triple Quadrupole Standards5
MS Single Quadrupole Standards7
MS Time-of-Flight Standards8
GC/MS Standards9
ACQUITY UPLC with UV Detection11
HPLC with UV Detection11
UPLC or HPLC with UV Detection Suitability Standards 12
HPLC with Fluorescence Detection
GC Standards
HPLC with RI, ELSD, Electrochemical,
or Conductivity Detection13
SFC and UPC ² Standards
pH Buffers14
Standards for Application Performance
Amino Acid Analysis17
Peptide Analysis
Glycan and Oligonucleotide Analysis
Protein Analysis
Small Molecule Analysis21
Environmental Analysis21
Food and Beverage Analysis
QC Trace Metal Standards
Ion Calibration Standards
Cation Analysis
Gel Permeation Chromatography (GPC) Standards 23
Gel Permeation Chromatography (GPC) Standards Kits 24
Gel Filtration Chromatography (GFC) Standards and Kits 24
Reagents
AccQ•Tag Ultra for UPLC Amino Acid Analysis
AccQ•Tag Amino Acid Analysis
Pico•Tag Amino Acid Analysis
Paired Ion Chromatography (PIC) Reagents
Food and Beverage Reagents
Reagents and Kits for Capillary Ion Analysis
J J J J J J J J J J J J J J J J J J J

Waters Standards and Reagents by Market Sector

Pharmaceutical	and Life	Sciences

Findi indicenticat and Life Sciences	
AccQ•Tag Ultra for UPLC and H-Class AAA Analysis	33
AccQ•Tag Amino Acid Analysis for HPLC AAA Analysis	34
Pico•Tag Amino Acid Analysis	34
Peptide Analysis	. 35
Protein Analysis	37
Glycan and Oligonucleotide Analysis	37
Small Molecule Analysis	. 38
Suitability Standards	. 38
Reagents and Kits for Capillary Ion Analysis	38
Paired Ion Chromatography (PIC) Reagents	39
Gel Filtration Chromatography (GFC)	
Standards and Kits	
pH Buffers	. 40
Food and Beverage	
Food and Beverage Analysis	. 43
Suitability Standards	. 43
Food and Beverage Reagents	. 43
pH Buffers	. 43
Environmental	
Reagents and Kits for Capillary Ion Analysis	47
SFC and UPC ² Standards	
Environmental Standards	47
Ion Calibration Standards	. 48
QC Trace Metal Standards	. 48
Cation Analysis	. 49
Suitability Standards	. 49
pH Buffers	. 49
Chemical Materials	
Gel Permeation Chromatography (GPC) Standards	53
Gel Permeation Chromatography (GPC) Standards Kits	53
SFC and UPC ² Standards	
Suitability Standards	. 54
pH Buffers	
Part Number Index	. 55
Alphabetical Index	. 57

Ordering Information61

System Performance Standards Key



Setup

Setup standards are used upon installation of a new system, start-up of a system which has been idle for a period of time, or to re-start a system that has been moved to a new location. These standards may include multiple components within a kit which allow step by step setup and calibration, or they can be single mixtures which will enable a full system assessment and initialization.



Calibration

These standards can be used for instrument calibration. Calibration is usually performed at a frequency that minimizes the risk of being out of calibration, in addition to major events such as when the instrument is moved and upon any major repair. Calibration can be confirmed as often as your standard operating procedures recommend. Look for solutions to calibrate your MS mass accuracy and resolution, UV detector linearity, UV and Fluorescence detector wavelength accuracy, Injector accuracy and precision, and solvent system flow rate accuracy.



Performance

Your system is comprised of many interdependent components working together. An issue with any one component can produce erroneous final results. All components performing correctly will produce results within an expected variability. Some of the system components include Mobile phase, Column, Tubing, Pumping, Injecting, Temperature controlling, Detecting, Data collection rate, and Integration parameters. An issue in any one of the components listed can affect your result.

Performance Check Standards should be utilized on a routine basis. This includes anytime a column is changed, a method is changed, the instrument is moved, or confirmation of system performance is desired. Confirming the performance of your system can be done on a system basis with suitability standards or per individual instrument component. Benchmark performance should be established after calibration and confirmed regularly.



Tuning

These Mass Spectrometry solutions are used routinely and are available for every ionization mode. Use the solution to optimize the mass spectrometer setting for optimal signal either manually or through Intellistart.





Standards for System Performance

Thousands of laboratories trust Waters to help them establish and document the performance of their instruments. We offer a broad range of products for most calibration needs, all backed by our reputation and professional support team.

- Calibrate and assess your systems' ability to perform your assays
- Perform a wide range of testing from set-up to calibration
- Monitor daily system proficiency to bring confidence to your assay execution
- Produce/generate inter/intra laboratory reproducibility

The following kits contain compounds appropriate to assess the performance or calibrate your triple quadrapole detector. Mass scale calibration, resolution, and sensitivity testing in different ionization modes can be performed using the compounds in these kits. They will help ensure optimal MS performance and bring confidence in your assay results.

MS Triple Quadrupole Standards

Calibration Performance Check Tuning

Description	Usage	Volume	Part No.
MS Set Up Solution, Atmospheric Pressure Ionization	S		700001594
2 ng/μL reserpine, 25 ng/μL PPG 1000, 50 ng/μL triacetyl-β-cyclodextrin in 50/50 acetonitrile/4 mM ammonium acetate		15 mL	
MS Calibration Solution, Atmospheric Pressure Ionization (API)	C		700001593
2 μg/μL sodium iodide, 50 ng/μL of cesium iodide in 50/50 isopropanol/water		15 mL	
MS Calibration Reference Standard—PEG 400	C		700004595
10 μg/mL PEG 400, 30.8 μg/mL ammonium acetate in 50/50 acetonitrile/water		10 mL	
SQD Set Up Solution	S		700003105
2.5 ng/μL PPG 1000, 5.0 ng/μL triacetyl-β-cyclodextrin, 100 pg/μL verapamil, 0.31 mg/mL ammonium acetat 1 ng/μL leucine enkaphalin in 50/50 acetonitrile/water	e,	10 mL	
Quadrapole: Tuning, Gradient, Sensitivity, and Precision Standard	SP T		186006861
200 pg/μL acetaminophen (4-acetamidophenol), 200 pg/μL caffeine, 100 pg/μL sulfadimethoxine, 50 pg/μL verapam 50 pg/μL chloramphenicol, 5000 pg/μL 17-α-hydroxyprogesterone in 90/10 methanol/water	il,	25 mL	
Reserpine Kit	P O		700001673
Reserpine Stock Solution: 100 pg/μL reserpine stock solution in 50/50 methanol/water Reserpine Working Solution: 1 pg/μL reserpine working solution in 50/50 methanol/water		10 mL 10 mL	
MS Solution—Daidzein (250 μg/L)	P		700002339
250 μg/L daidzein in 70/30 water/acetonitrile		10 mL	
Xevo® TQ-S Calibration and Resolution Standards	00		700005471
Calibration Solution: 0.1% ortho-phosphoric acid in 50/50 acetonitrile/water Intellistart™ Resolution Solution: 2 µg/µL sodium iodide, 50 ng/µL cesium iodide, 20 ng/µL methyl sulfate sodium salt in 50/50 2-propanol/water		30 mL	
Xevo TQ-S Qualification Standards Kit	G D	JO IIIL	700005472
Funing Solution 100X: 2 ng/μL acetominophen, 2 ng/μL caffeine, 1 ng/μL sulfadimethoxine, 500 pg/μL verapamil,			100003412
500 pg/μL chloramphenicol, 50 ng/μL 17-alpha-hydroxyprogesterone in 90/10 methanol/water Set-up Solution: 250 pg/μL PPG1000, 500 pg/μL triacetyl-β-cyclodextrin, 10 pg/μL verapmil, 100 pg/μL leucine en	kaphalin,	25 mL	
0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water IntelliStart Resolution Solution 20X: 2 μg/μL sodium iodide, 50 ng/μL cesium iodide,	,	30 mL	
20 ng/µL methyl sulfate, sodium salt in 50/50 2-propanol/water		30 mL	
Calibration Solution: 0.1% ortho-phosphoric acid in 50/50 acetonitrile/water on Sabre APCI Pos Sensitivity: 10 mg 17-alpha-hydroxyprogesterone		30 mL	
ESI Negative Sensitivity: 50 mg chloramphenicol		Solid Solid	
ESI Positive Sensitivity: 5 mg Reserpine		Solid	
Carryover: 5 mg sulfadimethoxine		Solid	
Scanwave Optimization: 0.1 mg [Glu1]-fibrinopeptide B		Solid	
Aqueous ammonia		3 mL	
Formic acid		1 mL	
(2) Formic Acid solvents		1 mL	
Kevo TQ-S Set Up Standard	00		700005331
Resolution Set-up and Calibration Standard: 250 pg/μL PPG1000, 500 pg/μL triacetyl-β-cyclodextrin, 10 pg/μL ven 100 pg/μL leucine enkaphalin, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water	rapmil,	30 mL	

MS Triple Quadrupole Standards (continued)

Description	Usage	Volume	Part No.
ACQUITY® MS Test Solution	P		700002728
(2) 1 mg/mL sulfadimethoxine in methanol		10 mL	
1 mg/mL terfenadine in acetonitrile		10 mL	
1 mg/mL reserpine in acetonitrile		10 mL	
1 mg/mL acetominophen in acetonitrile		10 mL	
1 mg/mL caffeine in acetonitrile		10 mL	
Formic acid		10 mL	
ACQUITY/Quattro micro™ or Quattro Premier™ MS Start Up Solution Kit	S		700002741
(2) 1 mg/mL sulfadimethoxine in methanol		10 mL	
1 mg/mL of each sulfadimethoxine, terfenadine, reserpine, acetominophen, and caffeine in acetonitrile		10 mL	
Formic acid		10 mL	
MS Qualification Standards Kit	S @ P		700002157
(2) 250 μg/L daidzein in 70/30 water/acetonitrile		10 mL	
Setup Solution: 2 ng/μL reserpine, 25 ng/μL PPG 1000 and 50 ng/μL triacetyl-β-cyclodextrin			
in 50/50 acetonitrile/4 mM ammonium acetate		15 mL	
Calibration Solution: 2 μg/μL sodium iodide and 50 ng/μL cesium iodide in 50/50 2-propanol/water		15 mL	
1 mg/L daidzein in 70/30 water/acetonitrile		10 mL	
Formic acid		10 mL	
Xevo TQ Standards Kit	S O P O		700004045
Setup Solution/Tuning Solution: 200 pg/µL acetominophen, 200 pg/µL caffeine, 100 pg/µL sulfadimethoxine,			
50 pg/µL verapamil, 50 pg/µL chloramphenicol, 5000 pg/µL 17- $lpha$ -hydroxyprogesterone in 90/10 methanol/wate	r	25 mL	
Set-up Solution QP: 2.5 ng/μL PPG 1000, 5.0 ng/μL triacetyl-β-cyclodextrin, 100 pg/μL verapmil,			
1.0 ng/µL leucine enkaphalinL, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water		30 mL	
Calibration Solution: 2 μg/μL sodium iodide, 50 ng/μL cesium iodide in 50/50 2-propanol/water		30 mL	
on Sabre APCI Pos Sensitivity: 10 mg 17-α-hydroxyprogesterone		Solid	
ESI Negative Sensitivity: 50 mg chloramphenicol		Solid	
ESI Positive Sensitivity: 5 mg reserpine		Solid	
Ion Sabre APCI Neg Sensitivity: 5 mg sulfadimethoxine		Solid	
Scanwave Optimization: 0.1 mg [Glu1]-fibrinopeptide B		Solid	
Aqueous ammonia		3 mL	
(2) Formic acid		1 mL	
Atmospheric Pressure Ionization Test Solution (ZMD, QTof I/II/API-US/Quattro micro, Platform II, Quattro II, LCT)	S @ P		700000889
Set up Solution: 2 ng/μL reserpine, 25 ng/μL PPG 1000, 50 ng/μL triacetyl-β- cyclodextrin in 50/50 acetonitrile/4 mM ammonium acetate		10	
		10 mL	
Calibration Solution: 2 µg/µL sodium iodide, 50 ng/µL cesium iodide, in 50/50 2-propanol/water		10 mL	
Horse Heart Myoglobin: 3.4 mg horse heart myoglobin in 4 mL amber screw top vial		Solid	
Raffinose: 5 ng/µL raffinose in 50/50 acetonitrile/water		10 mL	
Reserpine Stock Solution: 50 ng/μL reserpine stock solution in 50/50 methanol/water		10 mL	
D . W. I. C. I		10 mL	
p-nitrophenol Stock: 50 ng/μL p-nitrophenol stock solution in 50/50 methanol/water		10 mL	
p-nitrophenol Stock: 50 ng/µL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/µL p-nitrophenol working solution in 50/50 methanol/water	900	10 mL	700002646
p-nitrophenol Stock: 50 ng/µL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/µL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit	S 0 		700002646
p-nitrophenol Stock: 50 ng/μL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/μL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit SetUp Solution: 2.5 ng/μL PPG 1000, 5.0 ng/μL triacetyl-β-cyclodextrin, 100 pg/μL verapamil,	S @ P	10 mL	700002646
p-nitrophenol Stock: 50 ng/µL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/µL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit SetUp Solution: 2.5 ng/µL PPG 1000, 5.0 ng/µL triacetyl-β-cyclodextrin, 100 pg/µL verapamil, 1.0 ng/µL leucine enkaphalin, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water	S G P	10 mL	700002646
Reserpine Working Solution: 1 pg/μL reserpine working solution in 50/50 methanol/water p-nitrophenol Stock: 50 ng/μL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/μL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit SetUp Solution: 2.5 ng/μL PPG 1000, 5.0 ng/μL triacetyl-β-cyclodextrin, 100 pg/μL verapamil, 1.0 ng/μL leucine enkaphalin, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Resolution: 1.0 μg/μL PPG 2000 solution, 10.2 mg/mL ammonium acetate in 50/50 acetonitrile/water		10 mL 10 mL 10 mL	700002646
p-nitrophenol Stock: 50 ng/µL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/µL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit SetUp Solution: 2.5 ng/µL PPG 1000, 5.0 ng/µL triacetyl-β-cyclodextrin, 100 pg/µL verapamil, 1.0 ng/µL leucine enkaphalin, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Resolution: 1.0 µg/µL PPG 2000 solution, 10.2 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Mass Measurement Accuracy: 100.0 ng/µL PEG 1000, 0.17 mg/mL ammonium acetate in 50/50 acetonitrile/water		10 mL 10 mL 10 mL 10 mL	700002646
p-nitrophenol Stock: 50 ng/µL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/µL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit SetUp Solution: 2.5 ng/µL PPG 1000, 5.0 ng/µL triacetyl-β-cyclodextrin, 100 pg/µL verapamil, 1.0 ng/µL leucine enkaphalin, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Resolution: 1.0 μg/µL PPG 2000 solution, 10.2 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Mass Measurement Accuracy: 100.0 ng/µL PEG 1000, 0.17 mg/mL ammonium acetate in 50/50 acetonitrile/water SSI Sensitivity: 5 mg reserpine		10 mL 10 mL 10 mL 10 mL Solid	700002646
p-nitrophenol Stock: 50 ng/µL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/µL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit SetUp Solution: 2.5 ng/µL PPG 1000, 5.0 ng/µL triacetyl-β-cyclodextrin, 100 pg/µL verapamil, 1.0 ng/µL leucine enkaphalin, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Resolution: 1.0 μg/µL PPG 2000 solution, 10.2 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Mass Measurement Accuracy: 100.0 ng/µL PEG 1000, 0.17 mg/mL ammonium acetate in 50/50 acetonitrile/water SSI Sensitivity: 5 mg reserpine ESI Sensitivity: 10 mg raffinose		10 mL 10 mL 10 mL Solid Solid	700002646
p-nitrophenol Stock: 50 ng/µL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/µL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit SetUp Solution: 2.5 ng/µL PPG 1000, 5.0 ng/µL triacetyl-β-cyclodextrin, 100 pg/µL verapamil, 1.0 ng/µL leucine enkaphalin, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Resolution: 1.0 µg/µL PPG 2000 solution, 10.2 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Mass Measurement Accuracy: 100.0 ng/µL PEG 1000, 0.17 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Sensitivity: 5 mg reserpine ESI Sensitivity: 10 mg raffinose lon Sabre APCI Sensitivity: 10 mg 17-α-hydroxyprogesterone		10 mL 10 mL 10 mL Solid Solid Solid	700002646
p-nitrophenol Stock: 50 ng/μL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/μL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit SetUp Solution: 2.5 ng/μL PPG 1000, 5.0 ng/μL triacetyl-β-cyclodextrin, 100 pg/μL verapamil, 1.0 ng/μL leucine enkaphalin, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Resolution: 1.0 μg/μL PPG 2000 solution, 10.2 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Mass Measurement Accuracy: 100.0 ng/μL PEG 1000, 0.17 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Sensitivity: 5 mg reserpine ESI Sensitivity: 10 mg raffinose Ion Sabre APCI Sensitivity: 10 mg 17-α-hydroxyprogesterone ESI Sensitivity: 50 mg chloramphenicol		10 mL 10 mL 10 mL Solid Solid Solid Solid	700002646
p-nitrophenol Stock: 50 ng/μL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/μL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit SetUp Solution: 2.5 ng/μL PPG 1000, 5.0 ng/μL triacetyl-β-cyclodextrin, 100 pg/μL verapamil, 1.0 ng/μL leucine enkaphalin, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Resolution: 1.0 μg/μL PPG 2000 solution, 10.2 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Mass Measurement Accuracy: 100.0 ng/μL PEG 1000, 0.17 mg/mL ammonium acetate in 50/50 acetonitrile/wates ISI Sensitivity: 5 mg reserpine ESI Sensitivity: 10 mg raffinose Ion Sabre APCI Sensitivity: 10 mg 17-α-hydroxyprogesterone ESI Sensitivity: 50 mg chloramphenicol API Calibration Solution: 2 μg/μL sodium iodide, 50 ng/μL cesium iodide in 50/50 2-propanol/water		10 mL 10 mL 10 mL Solid Solid Solid Solid 10 mL	700002646
p-nitrophenol Stock: 50 ng/μL p-nitrophenol stock solution in 50/50 methanol/water p-nitrophenol Working Solutions: 2 pg/μL p-nitrophenol working solution in 50/50 methanol/water TQD/Quattro Premier/XE/Quattro micro Kit SetUp Solution: 2.5 ng/μL PPG 1000, 5.0 ng/μL triacetyl-β-cyclodextrin, 100 pg/μL verapamil, 1.0 ng/μL leucine enkaphalin, 0.31 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Resolution: 1.0 μg/μL PPG 2000 solution, 10.2 mg/mL ammonium acetate in 50/50 acetonitrile/water ESI Mass Measurement Accuracy: 100.0 ng/μL PEG 1000, 0.17 mg/mL ammonium acetate in 50/50 acetonitrile/wates ISI Sensitivity: 5 mg reserpine ESI Sensitivity: 10 mg raffinose Ion Sabre APCI Sensitivity: 10 mg 17-α-hydroxyprogesterone ESI Sensitivity: 50 mg chloramphenicol		10 mL 10 mL 10 mL Solid Solid Solid Solid	700002646

These multi-component kits and individual standards listed below are used for set-up, tuning, calibration, and sensitivity measurements on a MS Single Quadrapole detector.

MS Single Quadrupole Standards

Description	Usage	Volume	Part No.
MS Calibration Solution, Atmospheric Pressure Ionization (API)	•		700001593
2 μg/μL sodium iodide, 50 ng/μL of cesium iodide in 50/50 isopropanol/water		15 mL	
MS Calibration Reference Standard—PEG 400	•		700004595
10 µg/mL PEG 400, 30.8 µg/mL ammonium acetate in 50/50 acetonitrile/water		15 mL	
SQD Set Up Solution	S		700003105
2.5 ng/μL PPG 1000, 5.0 ng/μL triacetyl-β-cyclodextrin, 100 pg/μL verapamil, 0.31 mg/mL ammonium acetate, 1 ng/μL leucine enkaphalin in 50/50 acetonitrile/water		10 mL	
Quadrapole: Tuning, Gradient, Sensitivity, and Precision Standard	P 1		186006861
200 pg/μL acetaminophen (4-acetamidophenol), 200 pg/μL caffeine, 100 pg/μL sulfadimethoxine, 50 pg/μL verapamil, 50 pg/μL chloramphenicol, 5000 pg/μL 17-α-hydroxyprogesterone		25 mL	
Reserpine Kit	•		700001673
Reserpine Stock Solution: 100 pg/μL reserpine stock solution in 50/50 methanol/water Reserpine Working Solution: 1 pg/μL reserpine working solution in 50/50 methanol/water		10 mL 10 mL	
MS Solution—Daidzein (250 μg/L)	P		700002339
250 μg/L daidzein in 70/30 water/acetonitrile		10 mL	
Performance and Operation Solutions for MS	P		700001332
10 mg/L daidzein in 30/70 acetonitrile/water			
SQD Standards Kit	S G D		700003093
Set up Solution QP: 2.5 ng/μL PPG 1000, 5.0 ng/μL triacetyl-β-cyclodextrin, 100 pg/μL verapamil, 0.31 mg/mL ammonium acetate, 1 ng/μL leucine enkaphalin in 50/50 acetonitrile/water Stock Solution: 1 ng/μL reserpine, 2 ng/μL chloramphenicol in 50/50 acetonitrile/water Tuning Solution: 100 pg/μL reserpine, 200 pg/μL chloramphenicol in 50/50 acetonitrile/water (2) Working Solution: 10 pg/μL reserpine, 20 pg/μL chloramphenicol in 50/50 acetonitrile/water API Calibration Solution: 2 μg/μL sodium iodide, 50 ng/μL cesium iodide in 50/50 2-propanol/water Mobile Phase Blank: 50/50 acetonitrile/water		10 mL 10 mL 10 mL 10 mL 10 mL 10 mL	
MS Solution—Daidzein (1 mg/L)	P		700002340
1 mg/L daidzein in 70/30 water/acetonitrile		10 mL	
Atmospheric Pressure Photoionization Test Solution	P		700002752
1 mg/L acetaminophen in 70/30 water/methanol		5 mL	
Atmospheric Pressure Photoionization Sensitivity Standard	P		700002837
5 mg naphthalene in a 4 mL amber screw top vial		Solid	
IonSABRE™ Positive Ion Sensitivity Standard	P		700002857
5 mg 17-α-hydroxyprogesterone in a 4 mL amber screw top vial		Solid	
Cholesterol Sample for Atmospheric Pressure Photoionization Sensitivity Test	P		700002820
(2) 5 mg Cholesterol Standard in 4 mL amber screw top vials		Solid	
MS Formic Acid Solution	P		700002341
Formic Acid		10 mL	



These multi-component and individual compound kits are commonly used to perform a series of system-level performance tests for MS time-of-flight installation and as routine use standards and calibrants. ESI positive and negative sensitivity and precision, APCI sensitivity and precision can be assessed. Many protein and digestion calibration standards are also available.

MS Time-of-Flight Standards

Description Description	Usage	Volume	Part No.
Glu 1]-Fibrinopeptide B Standard	@ P		70000472
0.1 mg [Glu1]-fibrinopeptide-B in 4 mL amber vial			
1S Leucine Enkephalin Kit	@ P		70000245
? ng/µL leucine enkephalin in 50/50 acetonitrile/water % acetic acid in 50/50 acetonitrile/water Empty 12 mL amber vial with cap		5 mL 5 mL	
1S Reference Standard/Reverse Peptide Kit	G		70000508
mg Ser-Asp-Gly-Arg-Gly mg Gly-Arg-Gly-Asp-Ser		Solid Solid	
lass Scale and Resolution Set Up Solution	S P		70000391
.0 μg/μL sodium iodide, 0.05 μg/μL rubidium iodide in 50/50 2-propanol/water		10 mL	
1S Reference Standard/Horse Heart Myoglobin	P		70000474
2) 2.5 mg horse heart myoglobin		Solid	
MassPREP™ Enolase Digestion Standard	P		18600232
nmol tryptically digested yeast enolase		Solid	
MassPREP ADH Digestion Standard	P		18600232
nmol tryptically digested yeast alcohol dehydrogenase (ADH)		Solid	
Tof Standards Kit without Bovine	@ P		70000476
2) 3 mg leucine enkephalin 2 μg/μL sodium iodide in 50/50 2-propanol/water 2) 2.5 mg D(+)-raffinose 2) 0.1 mg [Glu1]-fibrinopeptide B 0.1 M sodium hydroxide 5) 10 mg α-cyano-4-hydroxycinnamic acid EG mix in 50/50 acetonitrile/0.2% formic acid EG MALDI 10 mg/mL 50/50 acetonitrile/water		Solid 25 mL 15 mL 4 mL 4 mL Solid 2 mL 2 mL	
Tof Qualification Standards Kit	@ P		70000327
2) leucine enkephalin, 3 mg 2 μg/μL sodium iodide in 50/50 2-propanol/ water Bovine insulin, 3 mg 2) 2.5 mg horse heart myoglobin 2) 2.5 mg D(+)-raffinose 2) 0.1 mg [Glu1]-fibrinopeptide B 0.1 M sodium hydroxide in water 2.52 mg hemoglobin (human) 5) α-cyano-4-hydroxycinnamic acid, 5 x 10 mg		Solid 25 mL Solid Solid Solid Solid Solid 4 mL Solid Solid	
EG Mix in 50/50 acetonitrile/0.2 % formic acid		2 mL	









TruView™ LCMS Certified Vials

TruView LCMS certified vials are for scientists working at low analyte concentrations, concerned with analyte loss due to adsorption, and require a clean vial for MS analytical results.

Find out more at www.waters.com/vials



MS Time-of-Flight Standards (continued)

Description	Usage	Volume	Part No.
LCT Premier/XE Standards Test Kit	P		700004027
0.1 M sodium hydroxide in water		10 mL	
10% formic acid in water		10 mL	
(2) 2-propanol/water 90/10		10 mL	
(5) 50/50 acetonitrile/water		20 mL	
Positive Ion Tests: 2 ng/µL leucine enkephalin in 50/50 acetonitrile/water		10 mL	
Positive Ion Tests: Leucine enkephalin, 2 mg		Solid	
Negative Ion Tests: 5 ng/µL raffinose in 50/50 acetonitrile/water		10 mL	
Negative Ion Tests: Raffinose, 2 mg		Solid	
Resolution Tests: Melittin, 0.3 mg		Solid	
(2) Resolution Tests: [Glu1]-fibrinopeptide, 0.1 mg		Solid	
Positive Ion APcI: 1 ng/ μ L 17- α -hydroxyprogesterone in 70/30 acetonitrile/water		10 mL	
Negative Ion APcI: 1 $ng/\mu L$ sulfadimethoxine in 70/30 acetonitrile/water		10 mL	
Positive Ion APPI: 1 ng/μL cholesterol in methanol		10 mL	
Chromatographic Signal to Noise: Reserpine, 2.5 mg		Solid	
Chromatographic Signal to Noise: Ammonium acetate, 5 g		Soild	
ACQUITY UPLC® + MS Test Mix for Precision: 2 ng/μL sulfadimethoxine in 10/90 acetonitrile/water		10 mL	
ACQUITY UPLC + MS Test Mix for Gradient Performance: 2 ng/μL each acetominophen, caffeine,			
sulfadimethoxine, reserpine, terfenadine, in 10/90 acetonitrile/water		10 mL	
Other supplies include:			
(2) 10 mL Class A volumetric flask & stoppers, (2) Pasteur pipettes, Latex pipette bulb			
MS ETD Reagent Kit	P		700005669
(5) 100 mg nitrosobenzene		Solid	
(5) 1,3-dicyanobenzene 100 mg		Solid	
MS ETD Standards Kit	P		700005666
0.5 mg substance-P acetate salt hydrate		Solid	
(4) LCMS Certified amber vials			

The GC/MS system can provide reliable results when calibrated and maintained using the standards in the kits below. The multi and individual components in the kits listed below are used for set up, tuning, calibration, and sensitivity measurements on your GC/MS detector.

GC/MS Standards

Description	Usage	Volume	Part No.
Resolution check on all GC/MS products	P		186006862
(2) Heptacosafluorotributylamine		1 mL	
CI Optimization on GCTs	P		186006863
(2) 2,4,6-tris(trifluoromethyl)-1,3,5-triazine		1 mL	
CI/FI Optimization and Tuning on GCTs	P		186006864
(2) Chloropentafluorobenzene		1 mL	
FI Optimization on GCTs	P		186006865
(3) Ampules of 700 µL perfluorotrimethylcyclohexane, 200 µL heptacosatrifluorobutylamine, 20 µL hexafluorobenzene, 20 µL pentafluoro-benzene, 10 µL chloropentafluorobenzene, 20 µL xylene, 20 µL acetone		990 µL	
Calibration/Tuning on Autospecs	0		186006866
Perfluorokerosene, 2 g			

GC/MS Standards (continued)

Description	Usage	Volume	Part No.
Quattro micro GC Standards Kit	@ P		700002877
Resolution: (2) Heptacosafluorotributylamine Calibration: (2) 2,4,6-tris(trifluoromethyl) EI MRM: (2) 10 pg/µL hexachlorobenzene DIP Sensitivity: (2) 2.5 ng/µL hexa-chlorobenzene CI Positive Ion Sensitivity: (2) 50 pg/µL benzophenone CI MRM: (2) 2 pg/µL benzophenone EI Positive Sensitivity Full Scan CI Negative Ion Sensitivity: (4) 1 pg/µL octafluoronaphthalene EI Sensitivity SIR: (2) 50 fg/µL octafluoronaphthalene DCI Probe Sensitivity: 2 ng/µL caffeine Other supplies: 20 x 32 mm glass screw top vials, 20 x 5.75" glass pasteur pipet, 2 mm latex pipet bulb		1 mL 1 mL 1 mL 1 mL 1 mL 1 mL 1 mL 1 mL	
GCT Premier (EI, CI, DCI, SIP) Standards Kit	0 0 0		700002878
Resolution: (2) Heptacosafluorotributylamine Tuning and Calibration: (2) 2,4,6-tris(trifluoromethyl)-1,3,5-triazine Tuning: (2) chloropentafluorobezene GC/MS Mass Measurement Accuracy: (2) 10 ng/µL methyl stearate in hexane DIP EI Positive Mass Measurement Accuracy: (2) 50 ng/µL methyl stearate in hexane GC-EI Positive Sensitivity: (2) 1 pg/µL hexachlorobenzene in hexane DRE Demonstration: (2) 2.5 ng/µL hexachlorobenzene in hexane GC-CI Positive Sensitivity: (2) 100 pg/µL benzophenone in hexane DIP CI Positive Mass Measurement Accuracy: (2) 200 ng/µL benzophenone in hexane GC-CI Negative Sensitivity: (2) 1 pg/µL octafluoronaphthalene in hexane DIP CI Negative Mass Measurement Accuracy: (2) 50 ng/µL octafluoronaphthalene in hexane DIP CI Negative Mass Measurement Accuracy: (2) 20 ng/µL caffeine in hexane Other supplies include: 20 x 32 mm glass screw top vials, 20 x 5.75" glass Pasteur pipettes, 2 mm latex pipette bulb		1 mL 1 mL 1 mL 1 mL 1 mL 1 mL 1 mL 1 mL	
GCT Premier (FI, FD) Standards Kit	000		700002881
Resolution: (2) Heptacosafluorotributylamine Tuning and Calibration: (2) chloropentafluorobenzene Tuning Solution: $3 \times 700 \mu L$ perfluorotrimethylcyclohexane, $200 \mu L$ heptacosa, $20 \mu L$ hexafluorobenzene, $20 \mu L$ pentafluoro-benzene, $10 \mu L$ chloropentafluorobenzene, $20 \mu L$ xylene, $20 \mu L$ acetone FI Positive Sensitivity: (2) $100 pg/\mu L$ hexadecane FD Positive Sensitivity: (2) $1 ng/\mu L$ cholesterol FI DIP Positive Mass Measurement Accuracy: (2) $50 \mu g/\mu L$ benzophenone Other supplies include: $20 \times 32 mm$ glass screw top vials, $20 \times 5.75 mass mass pasteur pipettes, 2 mm latex pipette bulb$		1 mL 990 µL 1 mL 1 mL 1 mL	
AutoSpec Premier Standards Test Kit	P		700004021
200 pg/μL methyl stearate in hexane 10 ng/μL methyl stearate in hexane 1.0 mg hexatriacontane 200 pg/μL anthraquinone in dichloromethane 100 fg/μL 2,3,7,8-TCDD, 1.0 fg/μL 2,3,7,8-(13C)-TCDD in nonane (2) Nonane Hexadecane Xylene Perfluorokerosene, 2 g Heptacosa (2) Hexane Other supplies include: 12 glass screw top vials, 12 x 5.75" glass pasteur pipettes, 2 mm latex pipette bulb		2 mL 2 mL Solid 2 mL 1 mL 2 mL 2 mL 2 mL Solid 1 mL	





S Setup Calibration Performance Check Tuning



The following standards and standards kits are used for the calibration and assessment of ACQUITY UPLC and HPLC system performance. Use of these standards and kits to benchmark and continuously monitor typical baseline performance will enhance productivity and increase the accuracy of results.

ACQUITY UPLC with UV Detection

Description	Usage	Volume	Part No.
ACQUITY UPLC/UV Qualification Standards Kit	G P		700002642
One each: 0.01 mg/mL, 0.03 mg/r 0.04 mg/mL, 0.05 mg/mL, 0.08 m 0.1 mg/mL, 0.16 mg/mL, 0.20 µg/	ng/mL,		
4.0 mg/mL caffeine in 90/10 water. 0.03 mg/mL thiourea, 0.110 mg/m	/acetonitrile	10 mL	
acetophenone in 90/10 water/aceto 4.0 µg/mL of each: 2-acetylfuran, a acetophenone, propiophenone, buty benzophenone, and valerophenone	onitrile acetanilide, ylparaben,	10 mL	
90/10 water/acetonitrile (2) 90/10 water/acetonitrile		10 mL 10 mL	
ACQUITY UPLC/UV Qualificatio Standards Kit for 3 or 5 mm Flo		TOTAL	700004380
One each: 0.03 mg/mL, 0.09 mg/r 0.12 mg/mL, 0.15 mg/mL, 0.24 m 0.3 mg/mL, 0.48 mg/mL, 0.6 µg/m 12.0 mg/mL caffeine in 90/10 water 12.0 µg/mL of each: 2-acetylfuran, acetophenone, propiophenone, buty benzophenone and valerophenone	ng/mL, nL, and n/acetonitrile acetanilide, glparaben,	10 mL	
90/10 water/acetonitrile (2) 90/10 water/acetonitrile		10 mL 10 mL	
ACQUITY UPLC/UV Start Up Standards Kits	6	TOTIL	700002669
0.10 mg/mL caffeine in 90/10 water 4.0 μg/mL each: 2-acetylfuran, ace acetophenone, propiophenone, buty benzophenone and valerophenone,	tanilide, Jlparaben,	10 mL	
90/10 water/acetonitrile		10 mL	
ACQUITY UPLC/UV Qualificatio Standards Kit 25 mm Flow Cell			700002846
One each: 0.004 mg/mL, 0.012 m 0.016 mg/mL, 0.02 mg/mL, 0.03 0.04 mg/mL, 0.064 mg/mL, 0.08 1.6 mg/mL caffeine in 90/10 water.	2 mg/mL, μg/mL, and	10 mL	
0.12 mg/mL thiourea, 0.044 mg/n acetophenone in 90/10 water/aceto 4.0 µg/mL of each: 2-acetylfuran, a acetophenone, propiophenone, buty	nL onitrile acetanilide,	10 mL	
benzophenone and valerophenone 90/10 water/acetonitrile (2) 90/10 water/acetonitrile		10 mL 10 mL	
ACQUITY UPLC/UV and ELSD Qualification Standards Kit	@ P		700002848
0.1600 mg/mL caffeine in 90/10 water/acetonitrile		10 mL	
AquaAnalysis System Start Up Solution	S		700003870

HPLC with UV Detection

Description	Usage	Volume	Part No.
Alliance® HT Qualification Standards Kit	G P		700002390
0.030 mg/mL uracil, 0.01 mg/mL ca 0.030 mg/mL uracil, 0.02 mg/mL ca 0.030 mg/mL uracil, 0.03 mg/mL ca 0.030 mg/mL uracil, 0.05 mg/mL ca	ffeine ffeine	10 mL 10 mL 10 mL 10 mL	
0.030 mg/mL uracil, 0.06 mg/mL ca 0.010 mg/mL uracil, 0.1 mg/mL methu 0.010 mg/mL uracil	ffeine	10 mL 10 mL 10 mL	
All standards in 30/70 methanol/wa			
Alliance with RI Detector Qualification Standards Test Kit	O P		700002391
•		101	
0.10 mg/mL caffeine 0.20 mg/mL caffeine		10 mL 10 mL	
0.50 mg/mL caffeine		10 mL	
(2) 1.00 mg/mL caffeine		10 mL	
0.10 mg/mL methyl paraben		10 mL	
19.1/80.9 methanol/water by		TOTTL	
weight, diluent only		10 mL	
4.00 mg/mL caffeine		10 mL	
All standards in 19.1/80.9 methand water by weight	ol/		
HPLC with UV Standards Kit	G P		700003741
0.030 mg/mL uracil, 0.010 mg/mL of 0.030 mg/mL uracil, 0.035 mg/mL of 0.030 mg/mL uracil, 0.06 mg/ml of 0.030 mg/mL uracil, 0.14 mg/mL ca 0.030 mg/mL uracil, 0.22 mg/mL ca 0.010 mg/mL uracil, 0.20 mg/mL caffeine 4.00 mg/mL caffeine 0.010 mg/mL of	caffeine caffeine L caffeine ffeine	10 mL 10 mL 10 mL 10 mL 10 mL 10 mL 10 mL 10 mL	
All standards in 70/30 water/methan	nol		
HPLC HT with UV Standards Kit	O P		700003742
0.030 mg/mL uracil, 0.010 mg/mL of 0.030 mg/mL uracil, 0.020 mg/mL of 0.030 mg/mL uracil, 0.030 mg/mL of 0.030 mg/mL uracil, 0.050 mg/mL of (3) 0.030 mg/mL uracil, 0.060 mg/mL (2) 0.010 mg/mL uracil 0.40 µg/mL caffeine 4.00 mg/mL caffeine One ampule of blank mobile phase	affeine affeine affeine	10 mL 10 mL 10 mL 10 mL 10 mL 10 mL 10 mL 10 mL	
All standards in 70/30 water/methan	nol		
Setup Calibration Per	formance	Check ①	Tuning

HPLC with **UV** Detection (continued)

Description	Usage	Volume	Part No.
HPLC UV Linearity Solutions Flow Cells for Analytical	@ P		700001076
0.030 mg/mL uracil, 0.010 mg/mL 0.030 mg/mL uracil, 0.035 mg/mL (2) 0.030 mg/mL uracil, 0.060 mg/m 0.030 mg/mL uracil, 0.110 mg/mL 0.030 mg/mL uracil, 0.140 mg/mL 0.030 mg/mL uracil, 0.220 mg/mL 0.010 mg/mL uracil, 0.10 mg/mL uracil, 0.10 mg/mL methyl paraben 0.030 mg/mL uracil, 0.170 mg/mL	caffeine nL caffeine caffeine caffeine caffeine	10 mL 10 mL 10 mL 10 mL 10 mL 10 mL 10 mL	
All standards in 30/70 methanol/wa	ater		
$\textbf{AutoPurification}^{\text{\tiny{M}}} \ \textbf{System Standa}$	rd 🕑		716000765
(3) 2500 μ g/mL thionin, 3000 μ g/m thioflavin, 2500 μ g/mL crystal viole 3/1 water/methanol		10 mL	
Gradient Valve Standard	P		WAT042876
0.56 g/L propyl paraben in methano	ol	5 mL	
UV Wavelength Accuracy (Erbium Perchlorate)	0		WAT042885
10.0 g/L erbium perchlorate in wate	r	10 mL	
UV Detector Linearity Standards	Kit 📵		WAT042881
One each: 5.0 mg/L, 10.0 mg/L, 15 20.0 mg/L, 25.0 mg/L, and 30.0 m propyl paraben in methanol	-	10 mL	
Two Cuvette Kit	O P		700004155
(2) Quartz cuvettes, 190 through 80	00 nm	Empty	
Six Cuvette Set	O P		700004156
(6) Quartz cuvettes, 190 through 80	00 nm	Empty	
Caffeine Standard	@ P		700003233
1.000 mg/mL caffeine in acetonitril	е	10 mL	
UV Detector Linearity Standards for 3 mm Flow Cell	Kit		WAT047634
One each: 15.0 mg/L, 30.0 mg/L, 460.0 mg/L, 75.0 mg/L, and 90.0 m paraben in methanol	_	10 mL	
UV Wavelength Accuracy Kit (Uracil)			WAT042879
15 mg/L uracil in methanol		10 mL	
UV Detector PQ Test Solution			WAT042887
10 mg/L acetophenone, 10 mg/L propiophenone, 10 mg/L butyropheno 3.2 g/L acetone in 60/40 methanol/		10 mL	

UPLC or HPLC with UV Detection – Suitability Standards

Description	Usage	Volume	Part No.
Neutral Suitability Standard	P		186006360
0.02 mg/L Thiourea, 0.25 mg/mL Naphthalene, 0.4 mg/mL Acenaphtha in 50/50 acetonitrile/water	alene	1 mL	
Reversed-Phase Suitability Standard	P		186006363
0.016 mg/mL uracil, 0.02 mg/mL butyl-paraben, 0.06 mg/mL naphthal 0.4 mg/mL propanol, 0.34 mg/mL di phthalate, 0.2 mg/mL acenaphthene, 0.1 mg/mL amitriptylene in 65/35 methanol/20 mM K ₂ HPO ₄ buffer pH	propyl	1 mL	
Preparative Chromatography Standard	P		186006703
5 mg/mL Diclofenac sodium salt, 5 mg/mL Diphenhydramine hydrochlo 5 mg/mL Flavone in DMSO	oride,	1 mL	
Dye Suitability Standard Kit	P		716000765
AutoPurification™ System Standard: (3) 2500 μg/mL thionin, 3000 μg/m thioflavin, 2500 μg/mL crystal violet 3/1 water/methanol		10 mL	

 $\label{eq:maintain_pour_fluorescence} \mbox{ Maintain your fluorescence detection HPLC system using the kits provided below.}$

HPLC with Fluorescence Detection

Description	Usage	Volume	Part No.
Alliance with Fluorescence Qualification Standards Test Kit	O P		700002753
(2) 0.5 pg/µL, 1.0 pg/µL, 5.0 pg/µL a in 80/20 acetonitrile/water One each:10.0 pg/µL, 2.5 ng/mL and 2.5 µg/mL, anthracene in 80/20 acetonitrile/water		10 mL	
Fluorescence Detector Standard Solution	P		700003694
5.0 pg/µL anthracene in 20/80 water/acetonitrile		1 mL	
Fluorescence Detector Performand Standard Solution	ce P		WAT047685
0.10 mg/L anthracene in 70/30 acetonitrile/water		10 mL	

Setup Calibration Performance Check Tuning



GC Standards

The gas chromatograph can be calibrated and assessed using the kit below. Kits are provided for use with FID and ECD.

GC

Description	Usage	Volume	Part No.
GC Performance Standard	P		700001892
1.0 mg/mL caffeine in methylene ch	loride	10 mL	
GC Qualification Standards Kit	P		700005094
(2) Isooctane (blank) One each: 1 μL/mL, 2 μL/mL, 4 μL/m 6 μL/mL, 8 μL/mL, 10 μL/mL 1-octar		10 mL	
in isooctane	,	10 mL	
1% dodecane in isooctane		10 mL	
0.035 ng/µL lindane		10 mL	
GC Headspace Qualification Standards Kit	P		700005095
(3) 0.1% o-xylene in acetone (6) 0.02% 1,2,4-trichlorobenzene		20 mL	
in acetone		20 mL	

The kits below can be used to provide calibration and assessment of various detector including Refractive Index, Conductivity, Evaporative Light Scattering, and Electrochemical detection.

HPLC with RI, ELSD, Electrochemical or Conductivity Detection

Description	Usage	Volume	Part No.
ELS Detector Standard	P		700002387
250 mg/L acetaminophen in 90/10 water/methanol		10 mL	
Refractive Index Performance and Operation Solution	@ P		700001713
2.0 g/L caffeine in 19/81 methano	ol/water	10 mL	
Electrochemical Detector PQ Solution	O P		700002286
1 mg/L acetaminophen in 90/10 water/methanol		10 mL	
Conductivity Detector OQ/PQ Solution	G P		700002309
Flow Cell Constant Solution: (2) 1 mM potassium chloride salt in w OQ/PQ Solution: 500 ppm sodium 1,000 ppm potassium bromide in w	chloride,	1 mL	



Ordering online has never been easier or more secure! Go to **www.waters.com/order** and you can:

- $\begin{tabular}{c} \hline \end{tabular} \begin{tabular}{c} Check out our latest discounts for online shopping. \end{tabular}$
- Set up wish lists for important, upcoming projects.
- Find out pricing and product availability quickly and easily.
- eMail your cart or wish list to other project or purchasing colleagues.



The kits below should be used to assess the performance of the SFC or ACQUITY UPC $^{2^{tot}}$ System. There are multi-component and individual component kits.

SFC and UPC² Standards

Description	Jsage	Volume	Part No.
Waters Prep 15/30 SFC System Test Mix and Internal Standard	P		700005675
21 mg antipyrin, 45 mg ketoprofen, 15 mg sulfamethazine, 60 mL container 125 mg acetazolamide, 125 mL contai	iner	Solid Solid	
Waters Prep 100 SFC System Test Mix and Internal Standard	P		700005674
6250 mg antipyrin, 5625 mg ketopro 625 mg fluorescein, 125 mL container 1250 mg sulfamethazine, 125 mL con		Solid Solid	
UPC ² Standard Mix	P		186006372
2 mg/mL each: 3-benzoylpyridine, cort 4-nitroaniline, 4,4'-biphenol in methar		1 mL	
UPC ² Gradient Standard	P		186006551
1 mg/mL coumarin, 1 mg/mL flavone, 2 mg/mL caffeine, 1 mg/mL thymine, 2 mg/mL prednisone in 2-propanol		1 mL	
UPC ² Caffeine Standard	P		186006614
1.0 mg/mL caffeine in 2-propanol		2 mL	
UPC ² Standard Ibuprofen	P		186006521
1.0 mg/mL ibuprofen in 2-propanol		2 mL	
UPC ² Standard Ketoprofen	P		186006522
1.0 mg/mL ketoprofen in 2-propanol		2 mL	
UPC ² Standard Flavone	P		186006523
1.0 mg/mL flavone in 2-propanol		2 mL	
UPC ² Standard Flurbiprofen	P		186006524
1.0 mg/mL flurbiprofen in 2-propanol		2 mL	

These pH Buffers are directly traceable to NIST SRMs, mercury free, guaranteed stable for at least one year after your receipt, and are supplied with a full certificate of analysis.

pH Buffers

Description	Volume	Part No.
pH 4 Liter		129
pH 4 Buffer	1 L	
pH 7 Liter		133
pH 7 Buffer	1 L	
pH 10 Liter		137
pH 10 Buffer	1 L	
pH 4 Pint		127
pH 4 Buffer	1 pint	
pH 7 Pint		131
pH 7 Buffer	1 pint	
pH 10 Pint		135
pH 10 Buffer	1 pint	

Performance Check



pH Buffers





Application Performance Standards & Reagents

Even the best laboratory equipment depends on accurate calibration. Waters understands this, and has developed a family of highly reliable, proven products for establishing and documenting instrument and application performance. Waters Systems and Application Performance Standards help ensure that your UPLC, GC, GC/MS, LC/MS, and SFC systems are capable of producing valid results day in and day out.

Many scientists use their instrumentation for very specific or targeted analyses. Waters has developed application performance standards to ensure that your systems will produce accurate and precise results according to your application needs.

- Confirmation of system set-up
- Confidence in methodology and data acquisition
- Streamline workflow by minimizing sample prep

Amino Acid Analysis

Description	Volume	Part No.
ACQUITY UPLC AAA Application Kit		176001279
This Kit is intended to enable existing ACQUITY UPLC Systems for AAA applications. Kit contains:		
Amino acid standard (4) sample tubes, 72/pk (3) Total recovery vials with caps Column stabilizer kit, 150 mm AccQ•Tag** Ultra derivitization kit AccQ•Tag Ultra C ₁₈ 1.7 µm, 2.1 x 100 mm Column AccQ•Tag Ultra Eluent A AccQ•Tag Ultra Eluent B Tube inlet .0025 ID PEEK** nut PDA assembly 2 µL Sample loop Column in-line filter kit UPLC AAA solution information set UPLC AAA app. solution startup tests Cert., AAA application & familiarization	1 L 1 L	
AccQ•Tag™ Ultra Chemistry Kit		176001235
This refill kit is intended to recharge the AccQ•Tag Ultra chemistries that are a part of the application add-on kit for ACQUITY UPLC Systems or ACQUITY H-Class Systems. This kit should not be purchased as part of an initial system. Kit includes: AccQ•Tag Ultra derivatization kit, 250 analyses AccQ•Tag Ultra C ₁₈ , 2.1 x 100 mm, 1.7 µm column AccQ•Tag Ultra Eluent A AccQ•Tag Ultra Eluent B (4) Sample tubes, 72/pk (10) Ampules of amino acid standard, hydrolysate (3) Total recovery vials, 100 vial/pk	950 mL 950 mL 1 mL	186003836 186003837 186003838 186003839 WAT007571 WAT088122 186000384C
AccQ•Tag Ultra Derivatization Kit		186003836
This reagent kit contains reagents to derivatize 250 samples and is a component of the UPLC Amino Acid Analysis Application Solution		
UPLC AAA H-Class Applications Kit		176002983
This kit is intended to enable existing ACQUITY UPLC H-Class systems for AAA applications. Kit includes:		
AccQ•Tag Ultra derivatization kit, 250 analyses AccQ•Tag Ultra C ₁₈ , 2.1 x 100 mm, 1.7 μm column AccQ•Tag Ultra Eluent A AccQ•Tag Ultra Eluent B (10) Ampules of amino acid standard, hydrolysate (3) Total recovery vials, 100 vials/pk Tube inlet 0.0025 ID PEEK nut PDA assembly Column in-line filter kit UPLC AAA H-Class solution information set AAA Application and Familiarization Service	950 mL 950 mL 1 mL	186003836 186003837 186003838 186003839 WAT088122 186000384C 430001783 205000343 716003230 741000299
AccQ●Tag Chemistry Kit		WAT052875
Reagent kit containing enough reagents for 250 analyses and includes:		
AccQ•Fluor Reagent 1, (5) vials AccQ•Fluor Reagent 2A, (5) vials AccQ•Fluor Reagent 2B, (5) vials AccQ•Fluor Reagent 3.9 x 150 mm (2) AccQ•Tag Eluent A, (4) Sample tubes, 72/pk (10) Ampules amino acid standard, hydrolysate	6 mL 3 mg 3 mL 1 L	WAT052885 WAT088122
Amino Acid Standard (AccQ•Tag, Pico•Tag®, AccQ•Tag Ultra)	7 1114	WAT088122
		W 11000122

The prepackaged life science standards eliminate the need to prepare, and store these materials for a convenient way to test and evaluate your LC or LC/MS applications.

Peptide Analysis

Description	Volume	Part No.
MassPREP Enolase Digestion Standard		186002325
Qualitative standard containing approximately 1 nmol of tryptically digested yeast enolase. Used for evaluation and benchmarking HPLC and UPLC Columns, evaluation of mass spectrometry and chromatographic instrument performance.	Solid	
MassPREP Phosphorylase b Standard		186002326
Qualitative standard containing approximately 1 nmol of tryptically digested phosphorylase b. Used for evaluation and benchmarking HPLC and UPLC Columns, evaluation of mass spectrometry and chromatographic instrument performance.	Solid	
MassPREP Bovine Hemoglobin Standard		186002327
Qualitative standard containing approximately 1 nmol of tryptically digested bovine hemoglobin. Used for evaluation and benchmarking HPLC and UPLC Columns, evaluation of mass spectrometry and chromatographic instrument performance.	Solid	
MassPREP ADH Digestion Standard		186002328
Qualitative standard containing approximately 1 nmol of tryptically digested yeast alcohol dehydrogenase (ADH). Used for evaluation and benchmarking HPLC and UPLC Columns, evaluation of mass spectrometry and chromatographic instrument performance.	Solid	
MassPREP BSA Digestion Standard		186002329
Qualitative standard containing approximately 1 nmol of tryptically digested bovine serum albumin (BSA). Used for evaluation and benchmarking HPLC and UPLC Columns, evaluation of mass spectrometry and chromatographic instrument performance.	Solid	
MassPREP Digestion Standard Kit		186002330
Kit contains one (1) each of the following tryptically digested standards:		
Yeast enolase Phosphorylase b Bovine hemoglobin Yeast alcohol dehydrogenase (ADH) Bovine serum albumin (BSA)	Solid Solid Solid Solid Solid	186002325 186002326 186002327 186002328 186002329
MassPREP Peptide Mixture		186002337
One vial with approximately 1 nmole of each:		
Allantoin (Vo Marker); RASG-1, angiotensin frag. 1-7, bradykinin; angiotensin II; angiotensin I, renin substrate, enolase T35, enolase T37, melittin. The peptide standard is useful to test LC columns and systems	درا: ۱	
dedicated to peptide separations.	Solid	

December 1	V-1	D N
Description	Volume	Part No.
MassPREP Peptide Mixture 5 Pack This is a 5 pack of 186002337. Each vial contains approximately 1 nmole of each:		186002338
Allantoin (Vo Marker); RASG-1, angiotensin frag. 1-7, bradykinin, angiotensin II, angiotensin I, renin substrate, enolase T35, enolase T37, melittin. The peptide standard is useful to test LC columns and systems dedicated to peptide separations.	Solid	
RapiGest™		186001860
The <i>Rapi</i> Gest SF (surfactant) radically enhances protein enzymatic digestions in terms of speed, higher quality, percent recovery and more reproducible sample preparation.		
Other formats available:		
RapiGest SF Surfactant, 1 mg/5 pack RapiGest SF Surfactant, 10 mg RapiGest SF Surfactant, 50 mg	Solid Solid Solid	186001861 186002123 186002122
Cytochrome C Digestion Standard		186006371
Qualitative digestion standard containing approximately 1 nmol of tryptically digested bovine heart cytochrome C, used for evaluation and benchmarking HPLC and UPLC Columns, evaluation of mass spectrometry and chromatographic instrument performance.	Solid	
Hi3 Phos B Standard		186006011
The Hi3 Phos B standard is a quantitative standard comprised of the top 6 ionizing peptides in the rabbit Phosphorylase B protein. This standard is primarily intended for use with the Hi3 quantification method for MSE proteomics data processed with ProteinLynx™ Global Server for samples of microbial origin. The standard may also be of use in the evaluation and benchmarking of proteomic LC/MS systems comprised of nanoACQUITY UPLC® and Synapt® and Xevo time-of-flight mass spectrometers.	Solid	
Hi3 E.Coli Standard		186006012
The Hi3 E.Coli standard is a quantitative standard comprised of the top 6 ionizing peptides in the E.coli ClpB protein. This standard is primarily intended for use with the Hi3 quantification method for MSE proteomics data processed with ProteinLynx Global Server for samples of animal origin. The standard may also be of use in the evaluation and benchmarking of proteomic LC/MS systems comprised of nanoACQUITY UPLC and Synapt and Xevo time-of-flight mass spectrometers.	Solid	

Peptide Analysis (continued)

Description	Volume	Part No.
MassPREP Digestion Standard Mix 1		186002865
Standard consists of a relatively close equimolar ratio of tryptically digested proteins:		
Yeast enolase Phosphorylase b Yeast alcohol dehydrogenase (ADH) Bovine serum albumin (BSA)	Solid	
Provides high confidence and coverage during expression profiling of complex protein samples		
MassPREP Digestion Standard Mix 2		186002866
Standard consists of a range of molar ratios for tryptically digested proteins:		
Yeast enolase Phosphorylase b Yeast alcohol dehydrogenase (ADH) Bovine serum albumin (BSA)	Solid	
Provides high confidence and coverage during expression profiling of complex protein samples		
MassPREP E. Coli Digest Standard		186003196
Qualitative standard containing approximately 1 nmol of tryptically digested and purified E.coli cytosolic protein.	Solid	
This standard provides a complex biological matrix for proteomics evaluations.		
MassPREP Phosphopeptide Standard Enolase		186003285
4 yeast enolase derived phosphorylated peptides: T18 1P, T19 1P, T43 1P, T43 2P.	Solid	
Used to optimize phosphopeptide detection in LC/MS, LC/UV, and MALDI-MS.		

Description	Volume	Part No.
MassPREP Enolase Digest with Phosphopeptides Mix		186003286
Yeast enolase spiked with 4 yeast enolase derived phosphorylated peptides: T18 1P, T19 1P, T43 1P, T43 2P.	Solid	
A more complex mixture used to optimize and troubleshoot phosphopeptide detection in LC/MS, LC/UV, and MALDI-MS.		
MassPREP Phosphopeptide Sample Kit—Enolase		186003287
(2) Vials:		
MassPREP enolase digestion standard MassPREP phosphopeptide	Solid	186002325
standard enolase	Solid	186003285
A kit that allows one to mix and optimize a complex standard per specific applications.		
MassPREP Enhancer (5 vials)		186003863
(5) 500 mg MassPREP Enhancer. A component in the MassPREP		
Phosphopeptide Enrichment Kit	Solid	
MassPREP Phosphopeptide Enrichment Kit		186003864
Phosphopeptide enrichment µElution plate		186003820
(2) MassPREP Enhancer MassPREP enolase digest	Solid	186003863
with phosphopeptides	Solid	186003286
This kit is developed fro selective enrichment of phosphopeptides from complex samples.		

Process Optimization

Reliable, precise qualitative and quantitative results in $10\ \text{minutes}.$

With Waters AccQ•Tag pre-column derivitization reagent, reversed-phase separation column, and tested methods, you can accurately separate and quantify amino acids derived from hydrolyzed proteins, cell culture and fermentation media, and other biological matrices.








Description	Volume	Part No.
MassPREP Glycoanalysis Kit		186002817
MassPREP HILIC µElution Plate (5) MassPREP MALDI Matrix DHB <i>Rapi</i> Gest SF	Solid Solid	186002780 186002333 186001860
Provides a simple and robust sample preparation without compromising sample recovery.		
MassPREP OST Standard		186004135
Approximately 1 nmol of 15, 20, 25, 30, and 35 nucleotide (nt) long oligodeoxythymidines, lyophilized. Useful in testing and confirming HPLC/UPLC, LC/MS, and column performance for oligo applications.	Solid	
Glycan Performance Test Standard		186006349
The Glycan Performance Test Standard is a 2-AB labeled human IgG-like standard that is QC verified to contain the components needed to benchmark and evaluate ACQUITY UPLC BEH Glycan, 1.7 µm Columns.	Solid	
Dextran Calibration Ladder		186006841
The Dextran Calibration Ladder is used to calibrate the HILIC column from retention time to GU values. This calibration ladder provides good peak shape and reliable identification		


2-AB labeled human IgG N-linked glycans using the Glycan Performance Test Standard

from 2 to 30 Glucose Units.

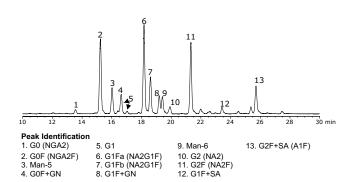

Solid

Typical chromatogram of 2-AB labeled human IgG N-linked glycans using the Glycan Performance Test Standard (Part No. 186003649)

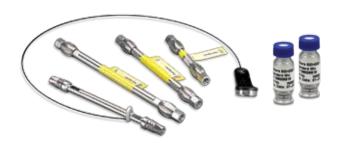

Protein Analysis

Description	Volume	Part No.
MassPREP Protein Standard Mix		186004900
Mixture of 6 intact proteins: enolase, BSA, ribonuclease A,cytochrome C, myoglobin and phos B, lyophilized. An attractive intact protein validation mixture that can be used for a variety of applications, specifically for benchmarking the Waters BEH C ₄ Column line.	Solid	
BEH125 SEC Protein Standard Mix		186006519
The BEH SEC Protein Standards are specifically designed to help aid in the benchmarking of the ACQUITY UPLC BEH125 SEC Column line. This mix contains 4 proteins: thyroglobulin, ovalbumin, ribonuclease A, and uracil.	Solid	
BEH200 SEC Protein Standard Mix		186006518
The BEH SEC Protein Standards are specifically designed to help aid in the benchmarking of the ACQUITY UPLC BEH200 SEC Column line. This mix contains 5 proteins: thyroglobulin, IgG, BSA, myoglobin, and uracil.	Solid	
BEH450 SEC Protein Standard Mix		186006842
The BEH SEC Protein Standards are specifically designed to help aid in the benchmarking of the ACQUITY UPLC BEH450 SEC Column line. This mix contains 5 proteins: thyroglobulin, IgG, BSA, myoglobin, and uracil.	Solid	
Intact mAb Mass Check Standard		186006552
The Intact mAb Mass Check Standard is an intact mouse protein purified by Protein-A with a known molecular weight. This standard is designed to be used for routine MS optimization, system setup and performance testing.	Solid	
IEX Anion Exchange Standard		186006869
For use with the Waters Protein-Pak™ Hi Res IEX Columns, this specific mix contains 3 proteins to allow a benchmark of the column when using anion exchange.	Solid	
IEX Cation Exchange Standard		186006870
For use with the Waters Protein-Pak Hi Res IEX Columns, this specific mix contains 3 proteins to allow a benchmark of the column when using cation exchange: bovine robonuclease A, cytochrome C, and chicken lysozyme.	Solid	

Small Molecule Analysis

Sillatt Plotecute Allatysis		
Description	Volume	Part No.
Analgesic Mix Standard		186006350
200 µg of each acetominophen, acetomidophenol, acetanilde, acetylsalicylic acid, caffeine, phenacetin, salicylic acid, lyophilized. Used for HPLC/UPLC method transfer and small molecule application work and validation.	Solid	
MetID Small Molecule Standard Mix		186006550
The MetID standard is a 3 component mix, omeprazole and 2 metabolites used for validation of the Waters metabonomic system and software.	Solid	
Neutral Suitability Standard		186006360
0.02 mg/L Thiourea, 0.25 mg/mL		
Naphthalene, 0.4 mg/mL Acenaphthalene in 50/50 acetonitrile/water	1 mL	
Reversed-Phase Suitability Standard		186006363
0.016 mg/mL uracil, 0.02 mg/mL butyl- paraben, 0.06 mg/mL naphthalene, 0.4 mg/mL propanol, 0.34 mgm/L dipropyl phthalate, 0.2 mg/mL acenaphthene, 0.1 mg/mL amitriptylene in 65/35 methanol/20mM K ₂ HPO ₄ buffer pH 7	1 mL	
Preparative Chromatography Standard		186006703
5 mg/mL Diclofenac sodium salt, 5 mg/mL Diphenhydramine hydrochloride, 5 mg/mL Flavone in DMSO	1 mL	
Dye Suitability Standard Kit		716000765
AutoPurification™ System Standard: (3) 2500 µg/mL thionin, 3000 µg/mL thioflavin, 2500 µg/mL crystal violet in 3/1 water/methanol	10 mL	



Environmental Analysis

Description	Volume	Part No.
20 Pesticide Mix Standard		186006348
10 mg/mL of each pesticide: atrazin, atrazin-desethyl, cyanazin, sebuthylazin, simazin, terbuthylazin, hexazinone, chlortoluron, diuron, linuron, metobromuron, metoxuron, monolinuron, metolachlor, methamidophos, acephate, methomyl, dicrotophos, buprofenzin, fenpropimorph in acetonitrile	1 mL	
Carbamate Analysis Standards Set		186004626
Carbamate/Carbamoxyloxine Pesticides, Potable Water: This 2 mL flame-sealed ampule yields up to 2 Liters after dilution. Each standard contains all analytes below in the required range 15-150 μg/L. Aldicarb, Aldicarb sulfone, Aldicarb sulfoxide, Baygon, Carbaryl, Carbfuran, 3-Hydroxycarbonfuran, Methiocarb, Methomyl, Oxamyl (vydate) Carbamate Pesticides, Priority Pollutant: This 2 mL flame-sealed ampule yields up to 2 Liters after dilution. Standard contains a random selection from the list below in the required range 5-200 μg/L. Aldicarb, Aldicarb sulfone, Aldicarb sulfoxide, Baygon, Carbaryl, Carbfuran, Difuron, 3-Hydroxycarbonfuran, Methiocarb, Methomyl, Oxamyl (vydate), Propham	2 mL	

Reliable SEC results depend on the quality of the reference materials used. Waters SEC calibration standards are precisely formulated to provide you with accurate molecular weight reference materials that are conveniently packaged to minimize errors in SEC calibration methods.

Find out more at www.waters.com/standards



Food and Beverage Analysis

Description	Volume	Part No.
Beverage Analysis Kit		176002534
Four each:		
Mobile phase Wash reagent Standards solution Solid standard	1 L 1 L 100 mL Solid	186006006 186006007 186006008 186006010
Beverage Analysis 5 Standards Solution		186006008
Caffeine 100 mg/L, benzoate 200 mg/L, sorbate 100 mg/L, acesulfame K 150 mg/L, saccharin 100 mg/L	100 mL	
Beverage Analysis Standard Solid		186006010
50 mg aspartame	Solid	
Soft Drink Analysis 5 Component Standard and Mobile Phase Kit		186003856
Soft Drink Analysis Mobile Phase: (4) 1 liter bottles containing a pre-mixed solution of denatured ethanol and water. Standards: 100 mg/L Caffeine, 200 mg/L benzoate, 100 mg/L sorbate, 150 mg/L acesulfame K in 7/93 ethanol/water	1 L	WAT039566
Aspartame (4) Aspartame, 100 mg	Solid	WAIU39566

Description	Volume	Part No.
Soft Drink Analysis 4 Component Standard		WAT036889
100 mg/L caffeine, 200 mg/L benzoate, 100 mg/L sorbate, 500 mg/L aspartame		
Beverage High Caffeine Standard		186006351
175 mg/L caffeine in 15/85 ethanol/water	500 mL	



Food and Beverage Analysis Kit

QC Trace Metal Standards

Description	Part No.
125 mL Aluminum	011
Aluminum in nitric acid, 1,000 mg/L	
125 mL Antimony	012
Antimony in nitric acid, 1,000 mg/L	
125 mL Arsenic	013
Arsenic in nitric acid, 1,000 mg/L	
125 mL Barium	014
Barium in nitric acid, 1,000 mg/L	
125 mL Beryllium	015
Beryllium in nitric acid, 1,000 mg/L	
125 mL Boron	016
Boron in nitric acid, 1,000 mg/L	
125 mL Cadmium	017
Cadmium in nitric acid, 1,000 mg/L	
125 mL Calcium	018
Calcium in nitric acid, 1,000 mg/L	
125 mL Chromium VI	019
Chromium VI in water, 1,000 mg/L	
125 mL Chromium	020
Chromium in nitric acid, 1,000 mg/L	

Description	Part No.
125 mL Cobalt	021
Cobalt in nitric acid, 1,000 mg/L	
125 mL Copper	022
Copper in nitric acid, 1,000 mg/L	
125 mL Iron	023
Iron in nitric acid, 1,000 mg/L	
125 mL Lead	024
Lead in nitric acid, 1,000 mg/L	
125 mL Magnesium	025
Magnesium in nitric acid, 1,000 mg/L	
125 mL Manganese	026
Manganese in nitric acid, 1,000 mg/L	
125 mL Mercury	027
Mercury in nitric acid, 1,000 mg/L	
125 mL Molybdenum	028
Molybdenum in nitric acid, 1,000 mg/L	
125 mL Nickel	029
Nickel in nitric acid, 1,000 mg/L	
125 mL Potassium	030
Potassium in nitric acid, 1,000 mg/L	

Description	Part No.
125 mL Selenium	031
Selenium in nitric acid, 1,000 mg/L	
125 mL Silicon	032
Silicon in nitric acid, 1,000 mg/L	
125 mL Silver	033
Silver in nitric acid, 1,000 mg/L	
125 mL Sodium	034
Sodium in nitric acid, 1,000 mg/L	
125 mL Strontium	035
Strontium in nitric acid, 1,000 mg/L	
125 mL Thallium	036
Thalium in nitric acid, 1,000 mg/L	
125 mL Tin	037
Tin in hydrochloric acid, 1,000 mg/L	
125 mL Titanium	038
Titanium in hydrochloric acid, 1,000 mg/L	
125 mL Vanadium	039
Vanadium in nitric acid, 1,000 mg/L	
125 mL Zinc	040
Zinc in nitric acid, 1,000 mg/L	

Ion Calibration Standards

Description	Part No.
125 mL Ammonia NH ₃	044
Ammonia as NH_3 in water, 1,000 mg/L	
125 mL Ammonia N	045
Ammonia as N in water, 1,000 mg/L	
125 mL Bromide	046
Bromide in water, 1,000 mg/L	
125 mL Chloride	047
Chloride in water, 1,000 mg/L	
125 mL Nitrate NO ₃	051-ERA
Nitrate as NO_3 in water, 1,000 mg/L	

Description	Part No.
125 mL Nitrate N	052-ERA
Nitrate as N in water, 1,000 mg/L	
125 mL Nitrite N	053
Nitrite as N in water, 1,000 mg/L	
125 mL Phosphate PO ₄	060
Phosphate as PO_4 in water, 1,000 mg/L	
125 mL Phosphate P	061-ERA
Phosphate as P in water, 1,000 mg/L	
125 mL Sulphate	062
Sulphate in water, 1,000 mg/L	

Description	Part No.
125 mL Phosphorus	063
Phosphorous in water, 1,000 mg/L	
125 mL Silica	064
Silica in water, 1,000 mg/L	
125 mL Bromate	065
Bromate in water, 1,000 mg/L	
125 mL Chlorate	066
Chlorate in water, 1,000 mg/L	
125 mL Chlorite	067
Chlorite in water, 1,000 mg/L	
125 mL Perchlorate	068
Perchlorate in water, 1,000 mg/L	

Cation Analysis

Description	Part No.
$125\ mL\ Ammonium\ NH_4$	78102
Ammonium as NH_4 in water, 100 mg/L	
125 mL Ammonium N	78104
Ammonium as N in water, 100 mg/L	
125 mL Lithium	78110
Lithium in water, 100 mg/L	

Description	Part No.
125 mL Calcium	K10
Calcium in water, 1,000 mg/L	
125 mL Magnesium	K11
Magnesium in water, 1,000 mg/L	
125 mL Potassium	K12
Potassium in water, 1,000 mg/L	

Description	Part No.
125 mL Sodium	K13
Sodium in water, 1,000 mg/L	
125 mL Nitrite NO ₂	K15
Nitrite as NO ₂ in water, 1,000 mg/L	

Gel Permeation Chromatography (GPC) Standards

Description	Part No.
Polystyrene Standard 400	WAT011590
Polystyrene, 400 Da, 10 g	
Polystyrene Standard 530	WAT011592
Polystyrene, 530 Da, 10 g	
Polystyrene Standard 950	WAT011594
Polystyrene, 950 Da, 10 g	
Polystyrene Standard 2,800	WAT011596
Polystyrene, 2,800 Da, 5 g	
Polystyrene Standard 6,400	WAT011598
Polystyrene, 6,400 Da, 5 g	
Polystyrene Standard 10,100	WAT011600
Polystyrene, 10,100 Da, 5 g	
Polystyrene Standard 17,000	WAT011602
Polystyrene, 17,000 Da, 5 g	
Polystyrene Standard 43,000	WAT011604
Polystyrene, 43,000 Da, 5 g	

Description	Part No.
Polystyrene Standard 110,000	WAT011606
Polystyrene, 110,000 Da, 5 g	
Polystyrene Standard 180,000	WAT011608
Polystyrene, 180,000 Da, 5 g	
Polystyrene Standard 430,000	WAT011612
Polystyrene, 430,000 Da, 5 g	
Polystyrene Standard 780,000	WAT011614
Polystyrene, 780,000 Da, 5 g	
Polystyrene Standard 1,300,000	WAT011616
Polystyrene, 1,300,000 Da, 1 g	
Polystyrene Standard 2,800,000	WAT011618
Polystyrene, 2,800,000 Da, 1 g	
Polystyrene Standard 3,600,000	WAT011620

Polystyrene, 3,600,000 Da, 1 g

Description	Part No.
Polystyrene Standard 4,300,000	WAT011622
Polystyrene, 4,300,000 Da, 1 g	
Polystyrene Standard 5,200,000	WAT011624
Polystyrene, 5,200,000 Da, 1 g	
Polystyrene Standard 6,200,000	WAT011626
Polystyrene, 6,200,000 Da, 1 g	
Polystyrene Standard 8,400,000	WAT011628
Polystyrene, 8,400,000 Da, 1 g	
Polystyrene Standard 20,000,000	WAT011630
Polystyrene, 20,000,000 Da, 1 g	


Gel Filtration Chromatography (GFC) Standards and Kits

For aqueous molecular weight determination.

ror aqueous motecutar weight determination.	
Description	Part No.
Polyethylene Oxide Standard 24,000	WAT011574
Polyethylene oxide, 24,000 Da, 0.5 g	
Polyethylene Oxide Standard 40,000	WAT011576
Polyethylene oxide, 40,000 Da, 0.5 g	
Polyethylene Oxide Standard 79,000	WAT011578
Polyethylene oxide, 79,000 Da, 0.5 g	
Polyethylene Oxide Standard 160,000	WAT011580
Polyethylene oxide, 160,000 Da, 0.5 g	
Polyethylene Oxide Standard 340,000	WAT011582
Polyethylene oxide, 340,000 Da, 0.5 g	
Polyethylene Oxide Standard 570,000	WAT011584
Polyethylene oxide, 570,000 Da, 0.5 g	
Polyethylene Oxide Standard 850,000	WAT011586
Polyethylene oxide, 850,000 Da, 0.5 g	
Polyethylene Glycol Standards Kit	WAT035711
1.0 g/vial polyethylene glycol at each molecular weight: 100, 200, 400, 600, 1,000, 1,500, 4,300, 7,000, 13,000, and 22,000	
Polyacrylic Acid Standards Kit	WAT035714
250 mg/vial polyacrylic acid at each molecular weight: 1,000, 3,000, 7,000, 15,000, 30,000, 70,000, 100,000, 300,000, 700,000, and 1,000,000	
Pullulan Kit	WAT034207
200 mg/vial Pullulan standards for aqueous MW testing: 5,000, 10,000, 20,000, 50,000, 100,000, 200,000, 400,000, and 800,000	
Polyethylene Oxide Kit	WAT011572
500 mg/vial Polyethylene oxide kit for aqueous MW testing: 24,000, 40,000, 79,000, 160,000, 340,000, 570,000, and 850,000	
Dextrans Standard	WAT054392
500 mg/vial dextrans at each molecular weight: 5,000, 12,000, 24,000, 48,000, 148,000, 273,000, 410,000, and 75,0000	

Gel Permeation Chromatography (GPC) Standards Kits

For non-aqueous molecular weight determination.	
Description	Part No.
Polymethylmethacrylate Mid MW Standards Kit	WAT035706
0.5 g/vial polymethylmethacrylate at each molecular weight: 2,400, 9,500, 31,000, 52,000, 10,0000, 170,000, 270,000, 490,000, 730,000, 1000000	
Polymethylmethacrylate Low MW Standards Kit	WAT035707
0.5 g/vial polymethylmethacrylate at each molecular weight: 1,000, 1,700, 2,500, 3,500, 5,000, 7,000, 10,000, 13,000, 20,000, 30,000	
Polyisoprene Standards Kit	WAT035708
0.5 g/vial polyisoprene at each molecular weight: 1000, 3,000, 10,000, 30,000, 70,000, 100,000, 300,000, 500,000, 1,000,000, 3,000,000	
Polybutadiene Standards Kit	WAT035709
0.5 g/vial polybutadiene at each molecular weight: 1,000, 3,000, 7,000, 10,000, 30,000, 70,000, 100,000, 300,000, 700,000, 1,000,000	
Polystyrene ReadyCal Standards 4 mL Kit	WAT058930
A complete kit of ready-to-use polystyrene calibration standards. Kit contains thirty 4 mL autosampler vials which contain four polystyrene standards per vial. There are three separate molecular weight ranges in each kit, ten units of each range. Range is from 400 to 2,000,000 Da.	
Polystyrene ReadyCal Standards 2 mL Kit	WAT058931
A complete kit of ready-to-use polystyrene calibration standards. Kit contains thirty 2 mL autosampler vials which contain four polystyrene standards per vial. There are three separate molecular weight ranges in each kit, ten units of each range. Range is from 400 to 2,000,000 Da.	
Polystyrene Kit Low-Mid MW	WAT011588
A complete kit of polystyrene at MW: 400, 530, and 950, 10 g each MW: 2,800, 6,400, 10,000, 17,000, 43,000, 110,000 and 180,000, 5 g each	
Polystyrene Kit Mid-High MW	WAT011610
A complete kit of polystyrene at MW: 430,000, 780,000, 5 g each MW: 1300,000, 2,800,000, 3,600,000, 4,300,000, 5,200,000, 6,200,000, 8,400,000, and 2,000,0000, 1 g each	
Polystyrene Kit SL-105	WAT034208
A complete kit of polystyrene at MW: 520, 950, 1,200, 1,800, 2,470, 3,770, 5,100, 7,600, 12,500, 17,000, 500 mg each	
Polystyrene Kit SM-105	WAT034209
A complete kit of polystyrene at MW: 1200, 3250, 10,200, 28,000, 68,000, 195,000, 490,000, 1,080,000, 1,750,000, and 2,750,000, 500 mg each	
Polystyrene Kit SH-75	WAT034210
A complete kit of polystyrene at MW: 450,000, 1,270,000, 2,300,000, 3,260,000, 4,340,000, 8,000,000, and 15,000,000, 500 mg each	

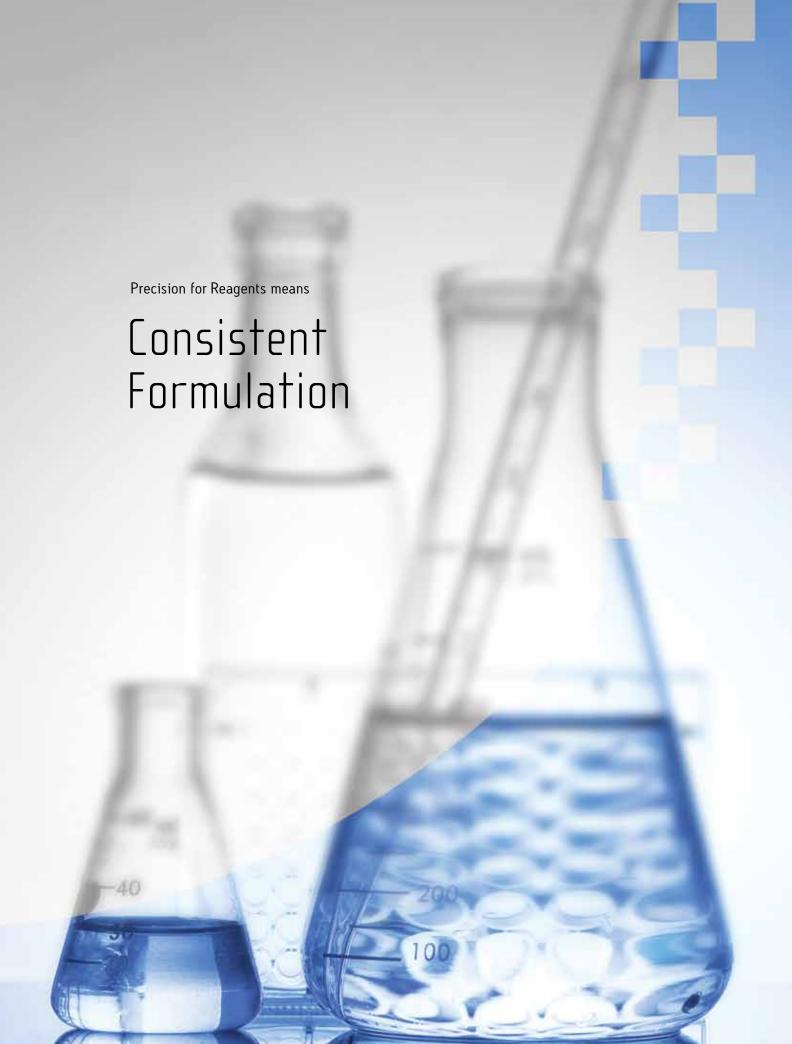



Reagents

The choice of mobile phases, solvents and additives can have a significant impact on chromatographic and spectrophotometric results. Waters provides standards kits and ready-to-use reagents that are prepared using the highest quality materials for all your critical analyses.

- Achieve optimum results from certified reagents
- Utilize precise formulation for ready-to-use capability
- Prepared using the highest quality materials

The AccQ•Tag Ultra chemistry products are a comprehensive, fully-tested set of reagents, columns, and eluents optimized for use with the UPLC Amino Acid Analysis Solution. This chemistry is based on Waters widely-used and well-understood AccQ•Tag derivatization method.

AccQ•Tag Ultra for UPLC and H-Class AAA Analysis

Description	Volume	Part No.
AccQ•Tag Ultra Chemistry Kit		176001235
The refill kit is intended to recharge the AccQ•Tag Ultra chemistries that are part of the application kit. This kit should be purchased for those that have already purchased the AccQ•Tag Ultra Application Solution. This kit is applicable to both ACQUITY UPLC and ACQUITY UPLC H-Class AAA Application Solutions. This kit should not be purchased as part of an initial system. Kit includes:		
AccQ•Tag Ultra derivatization kit, 250 analyses		186003836
AccQ•Tag Ultra C ₁₈ , 2.1 x 100 mm, 1.7 μm column	050 1	186003837
AccQ•Tag Ultra Eluent A concentrate	950 mL 950 mL	186003838 186003839
AccQ•Tag Ultra Eluent B (4) Sample tubes, 72/pk	950 IIIL	WAT007571
(10) Ampules of amino acid standard, hydrolysate	1 mL	WAT088122
(3) Total recovery vials, 100 vials/pk	I IIIL	186000384C
AccQ●Tag Ultra Eluent A		186003838
Concentrate	950 mL	
AccQ•Tag Ultra Eluent B		186003839
Concentrate	950 mL	
AccQ●Tag Ultra Derivatization Kit		186003836
Kit contains reagents to analyze 250 derivatized samples and is a component of the UPLC Amino Acid Analysis Application Solution.		

The HPLC AccQ•Tag method utilizes pre-column derivatized reagents that yield easily detected fluorescence adducts. The AccQ•Fluor™ reagent, 6-aminoquinolyl-N-hydroxysuccinimidyl carbamate derivatizes primary and secondary amines in a simple, single-step reaction to yield highly stable, fluorescent adducts. We offer the AccQ•Tag chemistry kits and individual reagents for your convenience.

AccQ•Tag Amino Acid Analysis for HPLC AAA Analysis

Description	Volume	Part No.
AccQ•Tag Chemistry Kit		WAT052875
Kit for up to 250 Analyses includes:		
AccQ•Fluor Reagent 1, (5) vials AccQ•Fluor Reagent 2A, (5) vials AccQ•Fluor Reagent 2B, (5) vials AccQ•Tag column, 3.9 x 150 mm	6 mL 3 mg 3 mL	
AccQ•Tag Eluent A concentrate (2)	1 L	
(4) Sample tubes, 72/pk (10) Ampules of amino acid standard, hydrolysate AccQ•Tag user guide	1 mL	WAT088122
Amino Acid Standard (AccQ•Tag, Pico•Tag, AccQ•Tag Ultra)		WAT088122
(10) Ampules of amino acid standards for use with AccQ•Tag, Pico•Tag, or AccQ•Tag Ultra Amino Acid Analysis Solutions	1 mL	
AccQ●Tag Eluent A		WAT052890
Concentrate	1 L	
AccQ●Tag Eluent B		WAT052895
Concentrate	11	
AccQ•Fluor Reagent Kit		WAT052880
Kit includes:		
AccQ•Fluor Reagent 1, (5) vials AccQ•Fluor Reagent 2A, (5) vials AccQ•Fluor Reagent 2B, (5) vials	6 mL 3 mg 4 mL	
The components of this kit are not available separately		

The Pico•Tag method provides a turnkey, guaranteed approach to modern HPLC amino acid analysis. Precolumn derivatization relies on the coupling reaction of the well known Edman Degradation. Purchase your individual reagents to recharge your analysis kit.

Pico•Tag Amino Acid Analysis

Description	Part No.
Pico•Tag Diluent	WAT088119
100 mL	
Pico•Tag Eluent 2	WAT010985
1 L	
Pico•Tag Eluent 2	WAT010965
(4) 1 L	
Pico•Tag Eluent 1	WAT010960
(4) 1 L	
Pico•Tag Eluent A	WAT088108
(4) 1 L	
Pico•Tag Eluent A for Low Pressure	186000689
(4) 1 L	
For sale in Japan only	
Pico•Tag Eluent B	WAT010983
1 L	
Pico•Tag Reagent Kit	WAT088123
Kit contains PITC, TEA, and standards	
Pico•Tag Chem Pkg Free AA	WAT091681
Designed for the analysis of physiologic amino acids.	
Free Amino Acid Chemistry Package includes:	
Free Amino Acid Analysis Column, 3.9 x 300 mm	WAT010950
Reagent kit	WATO10947 WATO10960
Pico•Tag Eluent A Pico•Tag Eluent 2	WATO 10965
Pico•Tag Diluent	WAT010303
Manual, column heater inserts, and sample tubes	
Pico•Tag Reagent Kit	WAT010947
Kit contains PITC, TEA, and standards A/N and B for Pico•Tag amino acid analysis of physiologic amino acids	
Pico•Tag Chemistry Package	WAT007360
Designed for the analysis of protein hydrolysates. Chemistry Package includes:	
Pico•Tag Column, 3.9 x 150 mm	WAT088131
Reagent kit (contains PITC, TEA, and Standards)	WAT088123
Pico•Tag Eluent A	WAT088108
Pico•Tag Eluent B Pico•Tag Diluent	WAT088112 WAT088119
Manual, column heater inserts, and sample tube	MUIOOOIIA
,	

Paired-ion chromatography is a technique for performing ion-exchange separations on a reversed-phase column. Hydrophobic paired-ion reagents added to the mobile phase are adsorbed onto the stationary phase. There they enhance the interaction of oppositely charged sample components with the stationary phase. Waters PIC Reagents are pre-mixed and buffered. Add 1 vial of reagent to 1 L of mobile phase and filter. Use the PIC reagents below for your paired ion chromatography technique.

Paired Ion Chromatography (PIC) Reagents

Description	Part No.
PIC A Low UV Reagent	WAT084189
5 vials of tetrabutylammonium hydrogen sulfate. Workable UV range in mobile phase is 200+ nm	
PIC B5 Low UV Reagent	WAT084198
5 vials of pentane sulfonic acid. Workable UV range in mobile phase is 200+ nm	
PIC B6 Low UV Reagent	WAT084199
5 vials of hexane sulfonic acid. Workable UV range in mobile phase is 200+ nm	
PIC B7 Low UV Reagent	WAT084282
5 vials of heptane sulfonic acid. Workable UV range in mobile phase is 200+ nm	
PIC B8 Low UV Reagent	WAT084283
5 vials of octane sulfonic acid. Workable UV range in mobile phase is 200+ nm	
PIC Reagent A	WAT085101
5 vials of tetrabutylammonium phosphate. Workable UV range in mobile phase is 240+ nm	
PIC Reagent B5	WAT085110
5 vials of pentane sulfonic acid. Workable UV range in mobile phase is 240+ nm	
PIC Reagent B6	WAT085140
5 vials of hexane sulfonic acid. Workable UV range in mobile phase is 240+ nm	
PIC Reagent B7	WAT085103
5 vials of heptane sulfonic acid. Workable UV range in mobile phase is 240+ nm	
PIC Reagent B8	WAT085142
5 vials of octane sulfonic acid. Workable UV range in mobile phase is 240+ nm	
PIC Reagent D4	WAT085466
5 vials of dibutylammonium phosphate. Workable UV range in mobile phase is 215+ nm	

Paired Ion Chromatography (PIC) Reagents (continued)

Description	Part No.
PIC Reagent Sample Kit (PIC A, B5, B6, B7, B8)	WAT085144
One each:	
PIC A	WAT085101
B5	WAT085110
B6	WAT085140
B7	WAT085103
B8	WAT085142
PITC Reagent	WAT088120
For use with Pico•Tag amino acid analysis method	
SAM I Reagent	WAT010873
(5) Reagent for carbohydrate analysis	

Waters application-specific reagents provide additional individual components for ongoing analysis.

Food and Beverage Reagents

Description	Volume	Part No.
Beverage Analysis Mobile-Phase Reagent		186006006
Acetate buffer with ethanol	1 L	
Beverage Analysis Wash Reagent		186006007
Ethanol solution	1 L	
Soft Drink Mobile Phase for SunFire™ Column Analysis		186003858
(4) Phosphate buffer with ethanol	1 L	
Soft Drink Mobile Phase for NovaPak® Column Analysis		WAT036888
(4) Potassium phosphate with ethanol	1 L	

The complete Beverage Analysis Kit can be found on page 41. The complete Soft Drink Kit can be found on page 41.

Waters reagents for capilliary ion analysis are electrolytes designed to make your ion analysis more convenient and reproducible. They also increase productivity in your lab by eliminating preparation time.

Reagents and Kits for Capillary Ion Analysis

Description	Volume	Part No.
CIA-Pak® MBT Solution		WAT036991
ElectroOsmotic flow modifier for working electrolyte		
IonSelect Anion Package		WAT049385
Waters IonSelect high mobility anion electrolyte Waters IonSelect high mobility anion test solution	1 L 5 mL	
IonSelect Cation Package		WAT049386
Waters IonSelect low mobility cation electrolyte Waters IonSelect low mobility cation test solution	1 L 5 mL	
CIA-Pak UV Cat 1 Reagent		WAT054780
Alpha hydroxyisobutyric acid 15 mL bottle for CIA system	100 mL	
IonSelect OFM Hydroxide Concentrate		WAT049387
IonSelect 100 mM OFM hydroxide concentrate		
CIA-Pak HIBA		WAT054775
2-hydroxyisobutyric acid, 2 g		
MS Formic Acid Solution		700002341
Formic Acid	10 mL	
Triethylamine		100000110
Triethylamine used in Pico•Tag Chem Pkg Free AA analyses	100 mL	

PIC Reagents are an easy solution to the problem of analyzing samples containing ionic species.

Learn more at www.waters.com/standards



MAKE YOUR ANALYSES, ROUTINE TESTINGS & VALIDATION EASIER

Save time with these 7 new standards. They are ready-to-use and specially formatted for their various applications.

Peptides



- Cytochrome C Digestion Standard (186006371): This standard is not only a valuable research standard but also a performance test standard.
- Hi3 Phos B Standard (186006011) and Hi3 E.coli Standard (186006012): The Hi3 Ecoli and Phos B standards are a specialized set of 6 synthetically prepared and highly purified peptides that provide the top 6 ionizing peptides of both the E.coli ClpB and rabbit Phosphorylase B proteins.

Proteins



- **ACQUITY UPLC® BEH200 SEC Protein Standard Mix (186006518):** The BEH200 SEC Protein Standard Mix provides a suitable mixture of proteins and an IgG to be able to provide a qualified test mix for the BEH200 SEC, 1.7 μm column line.
- ACQUITY UPLC BEH125 SEC Protein Standard Mix (186006519): The BEH125 Protein Standard Mix provides a suitable mixture of proteins within the molecular weight range of 1,000-80,000 daltons to provide a qualified text mix for the BEH125 SEC, 1.7 µm column line.
- Intact mAb Mass Check Standard (186006552): The Intact mAb Mass Check Standard is a LC/MS standard used as a qualitative tool for intact mass analysis and is designed for routine MS optimization, system setup and performance testing.

Glycans



■ Glycan Performance Test Standard (186006349): The Glycan Performance Standard is 2-AB labeled Human-like IgG standard provided to benchmark and measure life time performance of the ACQUITY UPLC BEH Glycan, 1.7 µm column line.

For more information visit www.waters.com/standards or contact your local sales specialist.



Waters





Pharmaceutical and Life Sciences

The measurement of therapeutic biological molecules requires standards and reagents designed specifically to ensure that your LC and LC/MS systems are fit for purpose. Waters has developed a portfolio of Life Science consumables designed for critical life science applications to drive successful results and add an additional layer of assurance.

- Choose from a range of specially designed products
- Achieve consistent, repeatable results
- Add confidence and assurance to your analyses
- Obtain optimum performance from your equipment by using the optimum standard.

The AccQ•Tag Ultra chemistry products are a comprehensive, fully-tested set of reagents, columns, and eluents optimized for use with the UPLC Amino Acid Analysis Solution. This chemistry is based on Waters widely-used and well-understood AccQ•Tag derivatization method.

AccQ•Tag Ultra for UPLC and H-Class AAA Analysis

Description	Volume	Part No.
ACQUITY UPLC AAA Application Kit		176001279
This Kit is intended to enable existing ACQUITY UPLC Systems for AAA applications. Kit contains:		
Amino acid standard (4) sample tubes, 72/pk (3) Total recovery vials with caps Column stabilizer kit, 150 mm AccQ•Tag" Ultra derivitization kit AccQ•Tag Ultra C ₁₈ 1.7 µm, 2.1 x 100 mm Column AccQ•Tag Ultra Eluent A AccQ•Tag Ultra Eluent B Tube inlet .0025 ID PEEK™ nut PDA assembly 2 µL Sample loop Column in-line filter kit UPLC AAA solution information set UPLC AAA app. solution startup tests Cert., AAA application & familiarization	1 L 1 L	
AccQ●Tag Ultra Chemistry Kit		176001235
The refill kit is intended to recharge the AccQ•Tag Ultra chemistries that are part of the application kit. This kit should be purchased for those that have already purchased the AccQ•Tag Ultra Application Solution. This kit is applicable to both ACQUITY UPLC and ACQUITY UPLC H-Class AAA Application Solutions. This kit should not be purchased as part of an initial system. Kit includes:		
AccQ•Tag Ultra derivatization kit, 250 analyses AccQ•Tag Ultra C ₁₈ , 2.1 x 100 mm, 1.7 μm column AccQ•Tag Ultra Eluent A concentrate AccQ•Tag Ultra Eluent B (4) Sample tubes, 72/pk (10) Ampules of amino acid standard, hydrolysate (3) Total recovery vials, 100 vials/pk	950 mL 950 mL 1 mL	186003836 186003837 186003838 186003839 WAT007571 WAT088122 186000384C
AccQ●Tag Ultra Eluent A		186003838
Concentrate	950 mL	
AccQ●Tag Ultra Eluent B		186003839
Concentrate	950 mL	
AccQ●Tag Ultra Derivatization Kit		186003836
Kit contains reagents to analyze 250 derivatized samples and is a component of the UPLC Amino Acid Analysis Application Solution.		
UPLC AAA H-Class Applications Kit		176002983
This kit is intended to enable existing ACQUITY UPLC H-Class systems for AAA applications. Kit includes:		
AccQ®Tag Ultra C18, 2.1 x 100 mm, 1.7 μm column AccQ®Tag Ultra Eluent A AccQ®Tag Ultra Eluent B (10) Ampules of amino acid standard, hydrolysate (3) Total recovery vials, 100 vials/pk Tube inlet 0.0025 ID PEEK nut PDA assembly Column in-line filter kit UPLC AAA H-Class solution information set AAA Application and Familiarization Service	950 mL 950 mL 1 mL	186003836 186003837 186003838 186003839 WAT088122 1860003840 430001783 205000343 716003230 741000299

The HPLC AccQ•Tag method utilizes pre-column derivatized reagents that yield easily detected fluorescence adducts. The AccQ•Fluor reagent, 6-aminoquinolyl-N-hydroxysuccinimidyl carbamate derivatizes primary and secondary amines in a simple, single-step reaction to yield highly stable, fluorescent adducts. We offer the AccQ•Tag chemistry kits and individual reagents for your convenience.

AccQ•Tag Amino Acid Analysis for HPLC AAA Analysis

Description	Volume	Part No.
AccQ●Tag Chemistry Kit		WAT05287
Kit for up to 250 Analyses includes:		
AccQ•Fluor Reagent 1, (5) vials AccQ•Fluor Reagent 2A, (5) vials AccQ•Fluor Reagent 2B, (5) vials AccQ•Tag column, 3.9 x 150 mm	6 mL 3 mg 3 mL	
AccQ•Tag Eluent A concentrate (2)	1 L	
(4) Sample tubes, 72/pk (10) Ampules of amino acid standard, hydrolysate AccQ•Tag user guide	1 mL	WAT08812
Amino Acid Standard (AccQ•Tag, Pico•Tag, AccQ•Tag Ultra)		WAT08812
(10) Ampules of amino acid standards for use with AccQ•Tag, Pico•Tag, or AccQ•Tag Ultra Amino Acid Analysis Solutions	1 mL	
AccQ•Tag Eluent A		WAT05289
Concentrate	1 L	
AccQ●Tag Eluent B		WAT05289
Concentrate	1 L	
AccQ•Fluor Reagent Kit		WAT05288
Kit includes:		
AccQ•Fluor Reagent 1, (5) vials AccQ•Fluor Reagent 2A, (5) vials AccQ•Fluor Reagent 2B, (5) vials	6 mL 3 mg 4 mL	
The components of this kit are not available separately		

The Pico Tag method provides a turnkey, guaranteed approach to modern HPLC amino acid analysis. Precolumn derivatization relies on the coupling reaction of the well known Edman Degradation. Purchase your individual reagents to recharge your analysis kit.

Pico•Tag Amino Acid Analysis

Description	Part No.
Pico•Tag Diluent	WAT088119
100 mL	
Pico•Tag Eluent 2	WAT010985
1 L	
Pico•Tag Eluent 2	WAT010965
(4) 1 L	
Pico•Tag Eluent 1	WAT010960
(4) 1 L	

Description	Part No.
Pico®Tag Eluent A	WAT088108
(4) 1 L	
Pico®Tag Eluent A for Low Pressure	186000689
(4) 1 L	
For sale in Japan only	
Pico®Tag Eluent B	WAT010983
1 L	
Pico•Tag Reagent Kit	WAT088123
Kit contains PITC, TEA, and standards	

Pico•Tag Amino Acid Analysis (continued)

Description	Part No.
Pico•Tag Chem Pkg Free AA	WAT091681
Designed for the analysis of physiologic amino acids.	
Free Amino Acid Chemistry Package includes:	
Free Amino Acid Analysis Column, 3.9 x 300 mm Reagent kit Pico•Tag Eluent A Pico•Tag Eluent 2 Pico•Tag Diluent Manual, column heater inserts, and sample tubes	WAT010950 WAT010947 WAT010960 WAT010965 WAT088119
Pico•Tag Reagent Kit	WAT010947
Kit contains PITC, TEA, and standards A/N and B for Pico•Tag amino acid analysis of physiologic amino acids	

Description	Part No.
Pico•Tag Chemistry Package	WAT007360
Designed for the analysis of protein hydrolysates. Chemistry Package includes:	
Pico•Tag Column, 3.9 x 150 mm Reagent kit (contains PITC, TEA, and Standards) Pico•Tag Eluent A Pico•Tag Eluent B Pico•Tag Diluent Manual, column heater inserts, and sample tube	WAT088131 WAT088123 WAT088108 WAT088112 WAT088119

The prepackaged life science standards eliminate the need to prepare, test, and store these materials for a convenient way to test and evaluate your LC or LC/MS applications.

Peptide Analysis

Volume	Part No.
	186002325
Solid	
	186002326
Solid	
	186002327
Solid	
Juliu	186002328
	100002320
Solid	
	Solid

Description	Volume	Part No.
MassPREP BSA Digestion Standard		186002329
Qualitative standard containing approximately 1 nmol of tryptically digested bovine serum albumin (BSA). Used for evaluation and benchmarking HPLC and UPLC Columns, evaluation of mass spectrometry and chromatographic instrument performance.	Solid	
MassPREP Digestion Standard Kit		186002330
Kit contains one (1) each of the following tryptically digested standards:		
Yeast enolase Phosphorylase b Bovine hemoglobin Yeast alcohol dehydrogenase (ADH) Bovine serum albumin (BSA)	Solid Solid Solid Solid Solid	186002325 186002326 186002327 186002328 186002329
MassPREP Peptide Mixture		186002337
One vial with approximately 1 nmole of each:		
Allantoin (Vo Marker); RASG-1, angiotensin frag. 1-7, bradykinin; angiotensin II; angiotensin I, renin substrate, enolase T35, enolase T37, melittin. The peptide standard is useful to test LC columns and systems dedicated to peptide separations.	Solid	
MassPREP Peptide Mixture 5 Pack		186002338
This is a 5 pack of 186002337. Each vial contains approximately 1 nmole of each:		
Allantoin (Vo Marker); RASG-1, angiotensin frag. 1-7, bradykinin, angiotensin II, angiotensin I, renin substrate, enolase T35, enolase T37, melittin. The peptide standard is useful to test LC columns and systems dedicated to peptide separations.	Solid	

Peptide Analysis (continued)

Description	Volume	Part No.
<i>Rapi</i> Gest [™]		186001860
The <i>Rapi</i> Gest SF (surfactant) radically enhances protein enzymatic digestions in terms of speed, higher quality, percent recovery and more reproducible sample preparation.		
Other formats available:		
<i>Rapi</i> Gest SF Surfactant, 1 mg/5 pack <i>Rapi</i> Gest SF Surfactant, 10 mg <i>Rapi</i> Gest SF Surfactant, 50 mg	Solid Solid Solid	186001861 186002123 186002122
Cytochrome C Digestion Standard		186006371
Qualitative digestion standard containing approximately 1 nmol of tryptically digested bovine heart cytochrome C, used for evaluation and benchmarking HPLC and UPLC Columns, evaluation of mass spectrometry and chromatographic instrument performance.	Solid	
Hi 3 Phos B Standard		186006011
The Hi3 Phos B standard is a quantitative standard comprised of the top 6 ionizing peptides in the rabbit Phosphorylase B protein. This standard is primarily intended for use with the Hi3 quantification method for MSE proteomics data processed with ProteinLynx™ Global Server for samples of microbial origin. The standard may also be of use in the evaluation and benchmarking of proteomic LC/MS systems comprised of nanoACQUITY UPLC® and Synapt® and Xevo time-of-flight		
mass spectrometers.	Solid	
Hi3 E.Coli Standard		186006012
The Hi3 E.Coli standard is a quantitative standard comprised of the top 6 ionizing peptides in the E.coli ClpB protein. This standard is primarily intended for use with the Hi3 quantification method for MSE proteomics data processed with ProteinLynx Global Server for samples of animal origin. The standard may also be of use in the evaluation and benchmarking of proteomic LC/MS systems comprised of nanoACQUITY UPLC and Synapt and Xevo time-of-flight mass spectrometers.	Solid	
MassPREP Digestion Standard Mix 1		186002865
Standard consists of a relatively close equimolar ratio of tryptically digested proteins:		
Yeast enolase Phosphorylase b Yeast alcohol dehydrogenase (ADH)	Solid	
Bovine serum albumin (BSA)		

Description	Volume	Part No.
MassPREP Digestion Standard Mix 2	- 120	186002866
Standard consists of a range of molar ratios for tryptically digested proteins:		
Yeast enolase Phosphorylase b Yeast alcohol dehydrogenase (ADH) Bovine serum albumin (BSA)	Solid	
Provides high confidence and coverage during expression profiling of complex protein samples		
MassPREP E. Coli Digest Standard		186003196
Qualitative standard containing approximately 1 nmol of tryptically digested and purified E.coli cytosolic protein.	Solid	
This standard provides a complex biological matrix for proteomics evaluations.		
MassPREP Phosphopeptide Standard Enolase		186003285
4 yeast enolase derived phosphorylated peptides: T18 1P, T19 1P, T43 1P, T43 2P.	Solid	
Used to optimize phosphopeptide detection in LC/MS, LC/UV, and MALDI-MS.		
MassPREP Enolase Digest with Phosphopeptides Mix		186003286
Yeast enolase spiked with 4 yeast enolase derived phosphorylated peptides: T18 1P, T19 1P, T43 1P, T43 2P.	Solid	
A more complex mixture used to optimize and troubleshoot phosphopeptide detection in LC/MS, LC/UV, and MALDI-MS.		
MassPREP Phosphopeptide Sample Kit—Enolase		186003287
(2) Vials:	0.1.1	10000005
MassPREP enolase digestion standard MassPREP phosphopeptide standard enolase	Solid Solid	186002325 186003285
A kit that allows one to mix and optimize a complex standard per specific applications.		
MassPREP Enhancer (5 vials)		186003863
(5) 500 mg MassPREP Enhancer. A component in the MassPREP Phosphopeptide Enrichment Kit	Solid	
MassPREP Phosphopeptide Enrichment Kit		186003864
Phosphopeptide enrichment µElution plate (2) MassPREP Enhancer MassPREP enolase digest with phosphopeptides	Solid Solid	186003820 186003863 186003286
This kit is developed fro selective enrichment of phosphopeptides from complex samples.		

Protein Analysis

Description	Volume	Part No.
MassPREP Protein Standard Mix		186004900
Mixture of 6 intact proteins: enolase, BSA, ribonuclease A,cytochrome C, myoglobin and phos B, lyophilized. An attractive intact protein validation mixture that can be used for a variety of applications, specifically for benchmarking the Waters BEH C ₄ Column line.	Solid	
BEH125 SEC Protein Standard Mix		186006519
The BEH SEC Protein Standards are specifically designed to help aid in the benchmarking of the ACQUITY UPLC BEH125 SEC Column line. This mix contains 4 proteins: thyroglobulin, ovalbumin, ribonuclease A, and uracil.	Solid	
BEH200 SEC Protein Standard Mix		186006518
The BEH SEC Protein Standards are specifically designed to help aid in the benchmarking of the ACQUITY UPLC BEH200 SEC Column line. This mix contains 5 proteins: thyroglobulin, IgG, BSA, myoglobin, and uracil.	Solid	
BEH450 SEC Protein Standard Mix	00114	186006842
The BEH SEC Protein Standards are specifically designed to help aid in the benchmarking of the ACQUITY UPLC BEH450 SEC Column line. This mix contains 5 proteins: thyroglobulin, IgG, BSA, myoglobin, and uracil.	Solid	10000042
Intact mAb Mass Check Standard		186006552
The Intact mAb Mass Check Standard is an intact mouse protein purified by Protein-A with a known molecular weight. This standard is designed to be used for routine MS optimization, system setup and performance testing.	Solid	
IEX Anion Exchange Standard		186006869
For use with the Waters Protein-Pak" Hi Res IEX Columns, this specific mix contains 3 proteins to allow a benchmark of the column when using anion exchange.	Solid	
IEX Cation Exchange Standard		186006870
For use with the Waters Protein-Pak Hi Res IEX Columns, this specific mix contains 3 proteins to allow a benchmark of the column when using cation exchange: bovine robonuclease A,		
cytochrome C, and chicken lysozyme.	Solid	

Performance-Qualified Rugged Methodologies that Assure Reproducible Results.

Find out more at www.waters.com/standards

Glycan and Oligonucleotide Analysis

D		D . N
Description	Volume	Part No.
MassPREP Glycoanalysis Kit		186002817
MassPREP HILIC µElution Plate		186002780
(5) MassPREP MALDI Matrix DHB	Solid	186002333
RapiGest SF	Solid	186001860
Provides a simple and robust sample preparation without compromising sample recovery.		
MassPREP OST Standard		186004135
Approximately 1 nmol of 15, 20, 25, 30, and 35 nucleotide (nt) long oligodeoxythymidines, lyophilized. Useful in testing and confirming HPLC/UPLC, LC/MS, and column performance for oligo applications.	Solid	
Glycan Performance Test Standard		186006349
The Glycan Performance Test Standard is a 2-AB labeled human IgG-like standard that is QC verified to contain the components needed to benchmark and evaluate ACQUITY UPLC		
BEH Glycan, 1.7 μm Columns.	Solid	
Dextran Calibration Ladder		186006841
The Dextran Calibration Ladder is used to calibrate the HILIC column from retention time to GU values. This calibration ladder provides good peak shape and reliable identification		
from 2 to 30 Glucose Units.	Solid	



Small Molecule Analysis

Description	Volume	Part No.
Analgesic Mix Standard		186006350
200 µg of each acetominophen, acetomidophenol, acetanilde, acetylsalicylic acid, caffeine, phenacetin, salicylic acid, lyophilized. Used for HPLC/UPLC method transfer and small molecule application work and validation.	Solid	
MetID Small Molecule Standard Mix		186006550
The MetID standard is a 3 component mix, omeprazole and 2 metabolites used for validation of the Waters metabonomic system and software.	Solid	

Suitability Standards

Description	Volume	Part No.
Neutral Suitability Standard		186006360
0.02 mg/L Thiourea, 0.25 mg/mL Naphthalene, 0.4 mg/mL Acenaphthalene in 50/50 acetonitrile/water	1 mL	
Reversed-Phase Suitability Standard		186006363
0.016 mg/mL uracil, 0.02 mg/mL butyl- paraben, 0.06 mg/mL naphthalene, 0.4 mg/mL propanol, 0.34 mgm/L dipropyl phthalate, 0.2 mg/mL acenaphthene, 0.1 mg/mL amitriptylene in 65/35 methanol/20mM K ₂ HPO ₄ buffer pH 7	1 mL	
Preparative Chromatography Standard		186006703
5 mg/mL Diclofenac sodium salt,		
5 mg/mL Diphenhydramine hydrochloride, 5 mg/mL Flavone in DMSO	1 mL	
Dye Suitability Standard Kit		716000765
AutoPurification™ System Standard: (3) 2500 μg/mL thionin, 3000 μg/mL thioflavin, 2500 μg/mL crystal violet in 3/1 water/methanol	10 mL	

Waters reagents for capilliary ion analysis are electrolytes designed to make your ion analysis more convenient and reproducible. They also increase productivity in your lab by eliminating preparation time.

Reagents and Kits for Capillary Ion Analysis

Description	Volume	Part No.
CIA-Pak® MBT Solution	Totume	WAT036991
		WAIO30991
ElectroOsmotic flow modifier for working electrolyte		
IonSelect Anion Package		WAT049385
Waters IonSelect high mobility anion electrolyte	1 L	
Waters IonSelect high mobility anion test solution	5 mL	
IonSelect Cation Package		WAT049386
Waters IonSelect low mobility cation electrolyte Waters IonSelect low mobility	1 L	
cation test solution	5 mL	
CIA-Pak UV Cat 1 Reagent		WAT054780
Alpha hydroxyisobutyric acid 15 mL bottle for CIA system	100 mL	
IonSelect OFM Hydroxide Concentrate		WAT049387
IonSelect 100 mM OFM hydroxide concentrate		
CIA-Pak HIBA		WAT054775
2-hydroxyisobutyric acid, 2 g		
MS Formic Acid Solution		700002341
Formic Acid	10 mL	
Triethylamine		100000110
Triethylamine used in Pico•Tag Chem Pkg Free AA analyses	100 mL	



Ordering online has never been easier or more secure! Go to **www.waters.com/order** and you can:

- Check out our latest discounts for online shopping.
 - Find out pricing and product availability quickly and easily.
- Set up wish lists for important, upcoming projects.
- eMail your cart or wish list to other project or purchasing colleagues.



Paired-ion chromatography is a technique for performing ion- exchange separations on a reversed-phase column. Hydrophobic paired-ion reagents added to the mobile phase are adsorbed onto the stationary phase. There they enhance the interaction of oppositely charged sample components with the stationary phase. Waters PIC Reagents are pre-mixed and buffered. Add 1 vial of reagent to 1L of mobile phase and filter. Use the PIC reagents below for your paired ion chromatography technique.

Paired Ion Chromatography (PIC) Reagents

3 1 3 ()	•
Description	Part No.
PIC A Low UV Reagent	WAT084189
5 vials of tetrabutylammonium hydrogen sulfate. Workable UV range in mobile phase is 200+ nm	
PIC B5 Low UV Reagent	WAT084198
5 vials of pentane sulfonic acid. Workable UV range in mobile phase is 200+ nm	
PIC B6 Low UV Reagent	WAT084199
5 vials of hexane sulfonic acid. Workable UV range in mobile phase is 200+ nm	
PIC B7 Low UV Reagent	WAT084282
5 vials of heptane sulfonic acid. Workable UV range in mobile phase is 200+ nm	
PIC B8 Low UV Reagent	WAT084283
5 vials of octane sulfonic acid. Workable UV range in mobile phase is 200+ nm	
PIC Reagent A	WAT085101
5 vials of tetrabutylammonium phosphate. Workable UV range in mobile phase is 240+ nm	
PIC Reagent B5	WAT085110
5 vials of pentane sulfonic acid. Workable UV range in mobile phase is 240+ nm	
PIC Reagent B6	WAT085140
5 vials of hexane sulfonic acid. Workable UV range in mobile phase is 240+ nm	

Description	Part No.
PIC Reagent B7	WAT085103
5 vials of heptane sulfonic acid. Workable UV range in mobile phase is 240+ nm	
PIC Reagent B8	WAT085142
5 vials of octane sulfonic acid. Workable UV range in mobile phase is 240+ nm	
PIC Reagent D4	WAT085466
5 vials of dibutylammonium phosphate. Workable UV range in mobile phase is 215+ nm	
PIC Reagent Sample Kit (PIC A, B5, B6, B7, B8)	WAT085144
One each:	
PIC A B5 B6 B7 B8	WAT085101 WAT085110 WAT085140 WAT085103 WAT085142
PITC Reagent	WAT088120
For use with Pico•Tag amino acid analysis method	
SAM I Reagent	WAT010873
(5) Reagent for carbohydrate analysis	

ACQUITY UPLC BEH125, BEH200 and BEH450 SEC, columns and standards are ideally suited for the separation and analysis of peptides and proteins including monoclonal antibody monomers from aggregates.

Learn more at www.waters.com/biosep



Gel Filtration Chromatography (GFC) Standards and Kits

For aqueous molecular weight determination.

•	
Description	Part No.
Polyethylene Oxide Standard 24,000	WAT011574
Polyethylene oxide, 24,000 Da, 0.5 g	
Polyethylene Oxide Standard 40,000	WAT011576
Polyethylene oxide, 40,000 Da, 0.5 g	
Polyethylene Oxide Standard 79,000	WAT011578
Polyethylene oxide, 79,000 Da, 0.5 g	
Polyethylene Oxide Standard 160,000	WAT011580
Polyethylene oxide, 160,000 Da, 0.5 g	
Polyethylene Oxide Standard 340,000	WAT011582
Polyethylene oxide, 340,000 Da, 0.5 g	
Polyethylene Oxide Standard 570,000	WAT011584
Polyethylene oxide, 570,000 Da, 0.5 g	
Polyethylene Oxide Standard 850,000	WAT011586
Polyethylene oxide, 850,000 Da, 0.5 g	
Polyethylene Glycol Standards Kit	WAT035711
1.0 g/vial polyethylene glycol at each molecular weight: 100, 200, 400, 600, 1,000, 1,500, 4,300, 7,000, 13,000, and 22,000	
Polyacrylic Acid Standards Kit	WAT035714
250 mg/vial polyacrylic acid at each molecular weight: 1,000, 3,000, 7,000, 15,000, 30,000, 70,000, 100,000, 300,000, 700,000, and 1,000,000	
Pullulan Kit	WAT034207
200 mg/vial Pullulan standards for aqueous MW testing: 5,000, 10,000, 20,000, 50,000, 100,000, 200,000, 400,000, and 800,000	
Polyethylene Oxide Kit	WAT011572
500 mg/vial Polyethylene oxide kit for aqueous MW testing: 24,000, 40,000, 79,000, 160,000, 340,000, 570,000, and 850,000	
Dextrans Standard	WAT054392
500 mg/vial dextrans at each molecular weight: 5,000, 12,000, 24,000, 48,000, 148,000, 27,000, 41,000,	

These pH Buffers are directly traceable to NIST SRMs, mercury free, guaranteed stable for at least one year after your receipt, and are supplied with a full certificate of analysis.


pH Buffers

Description	Volume	Part No.
pH 4 Liter		129
pH 4 Buffer	1 L	
pH 7 Liter		133
pH 7 Buffer	1 L	
pH 10 Liter		137
pH 10 Buffer	1 L	
pH 4 Pint		127
pH 4 Buffer	1 pint	
pH 7 Pint		131
pH 7 Buffer	1 pint	
pH 10 Pint		135
pH 10 Buffer	1 pint	



pH Buffers

273,000, 410,000, and 75,0000

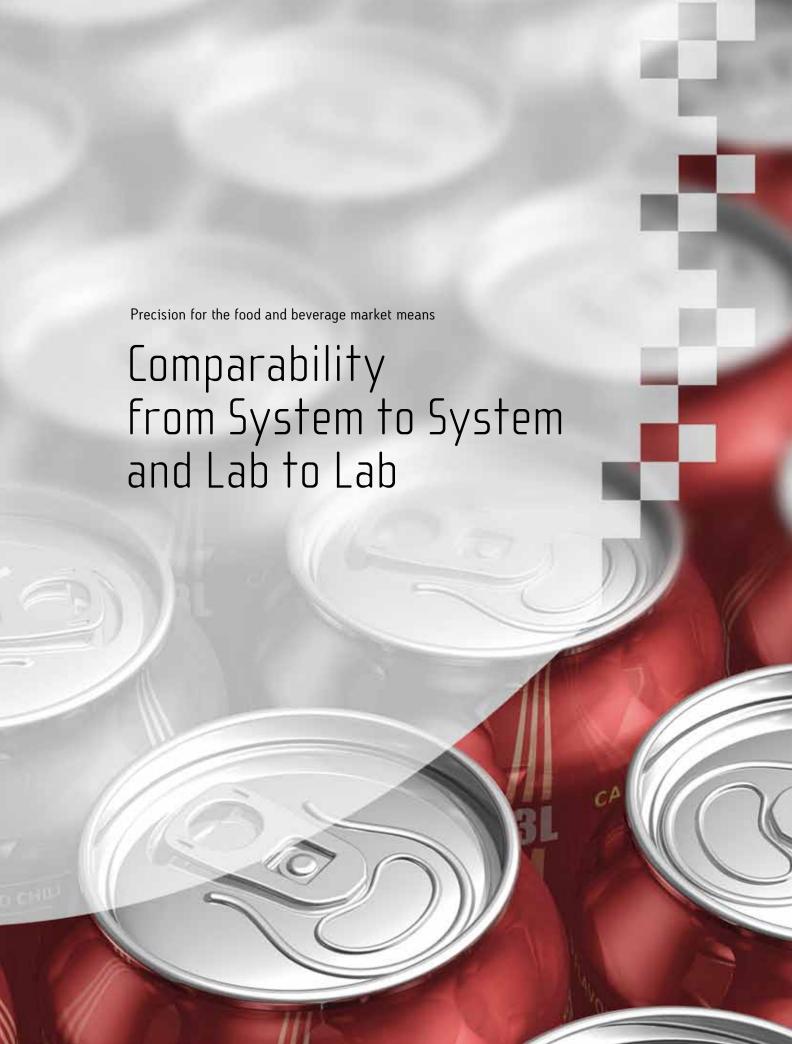



Food and Beverage

From QC testing of final product, to meeting the demand for new products with additional nutritional value, food manufacturers face difficult challenges. Waters can help meet those challenges with standards and reagents that offer the sensitivity, reproducibility and versatility required to ensure you can profitably manage laboratory requirements.

- Test complex food and beverage samples consistently and routinely
- Optimize analytical capability by improving confidence in results and instrumentation
- Take advantage of Waters' unique commitment and experience in the food industry

Food and Beverage Analysis

Description	Volume	Part No.
Beverage Analysis Kit		176002534
Four each:		
Mobile phase Wash reagent Standards solution Solid standard	1 L 1 L 100 mL Solid	186006006 186006007 186006008 186006010
Beverage Analysis 5 Standards Solution		186006008
Caffeine 100 mg/L, benzoate 200 mg/L, sorbate 100 mg/L, acesulfame K 150 mg/L, saccharin 100 mg/L	100 mL	
Beverage Analysis Standard Solid		186006010
50 mg aspartame	Solid	
Soft Drink Analysis 5 Component Standard and Mobile Phase Kit		186003856
Soft Drink Analysis Mobile Phase: (4) 1 liter bottles containing a pre-mixed solution of denatured ethanol and water. Standards: 100 mg/L Caffeine, 200 mg/L benzoate, 100 mg/L sorbate, 150 mg/L acesulfame K in 7/93 ethanol/water	1 L	
Aspartame (4) Aspartame, 100 mg	Solid	WAT039566
Soft Drink Analysis 4 Component Standard		WAT036889
100 mg/L caffeine, 200 mg/L benzoate, 100 mg/L sorbate, 500 mg/L aspartame		
Beverage High Caffeine Standard		186006351
175 mg/L caffeine in 15/85 ethanol/water	500 mL	

These pH Buffers are directly traceable to NIST SRMs, mercury free, guaranteed stable for at least one year after your receipt, and are supplied with a full certificate of analysis.

pH Buffers

Description	Volume	Part No.
pH 4 Liter		129
pH 4 Buffer	1 L	
pH 7 Liter		133
pH 7 Buffer	1 L	
pH 10 Liter		137
pH 10 Buffer	1 L	
pH 4 Pint		127
pH 4 Buffer	1 pint	
pH 7 Pint		131
pH 7 Buffer	1 pint	
pH 10 Pint		135
pH 10 Buffer	1 pint	

Suitability Standards

•		
Description	Volume	Part No.
Neutral Suitability Standard		186006360
0.02 mg/L Thiourea, 0.25 mg/mL Naphthalene, 0.4 mg/mL Acenaphthalene in 50/50 acetonitrile/water	1 mL	
Reversed-Phase Suitability Standard		186006363
0.016 mg/mL uracil, 0.02 mg/mL butyl- paraben, 0.06 mg/mL naphthalene, 0.4 mg/mL propanol, 0.34 mgm/L dipropyl phthalate, 0.2 mg/mL acenaphthene, 0.1 mg/mL amitriptylene in 65/35 methanol/20mM K ₂ HPO ₄ buffer pH 7	1 mL	
Preparative Chromatography Standard		186006703
5 mg/mL Diclofenac sodium salt, 5 mg/mL Diphenhydramine hydrochloride, 5 mg/mL Flavone in DMSO	1 mL	
Dye Suitability Standard Kit		716000765
AutoPurification™ System Standard: (3) 2500 μg/mL thionin, 3000 μg/mL thioflavin, 2500 μg/mL crystal violet in 3/1 water/methanol	10 mL	

Food and Beverage Reagents

Description	Volume	Part No.
Beverage Analysis Mobile-Phase Reagent		186006006
Acetate buffer with ethanol	1 L	
Beverage Analysis Wash Reagent		186006007
Ethanol solution	1 L	
Soft Drink Mobile Phase for SunFire™ Column Analysis		186003858
(4) Phosphate buffer with ethanol	1 L	
Soft Drink Mobile Phase for NovaPak® Column Analysis		WAT036888
(4) Potassium phosphate with ethanol	1 L	



Waters Beverage Analysis Kit is a comprehensive solution for determination of soft drink components. Available for use with UPLC and HPLC.

Learn more at www.waters.com/standards


We Are Raising Your Standards



Waters understands that the quality of the standards and reagents you use directly correlates to the quality of your results and final products. That is why we have made analytical standards and reagents our business. Our seal symbolizes Waters commitment to you.

waters.com/standards



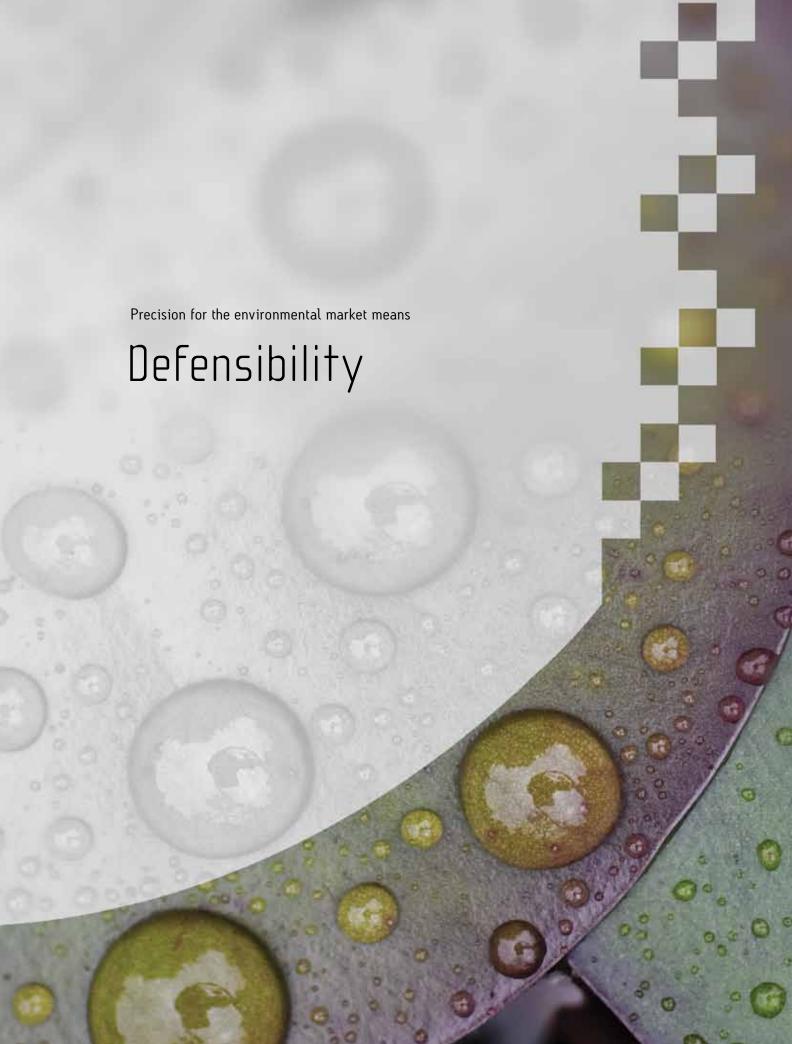



Environmental

Environmental quality issues are complex, challenging, and ever expanding. Understanding the complex nature of this type of analysis drove Waters' commitment to partner for a cleaner environment. Waters instrumentation, standards and sample preparation products are used for environmental analysis in laboratories throughout the world.

- Help ensure safety and compliance
- Increase laboratory efficiency
- Improve quality of results

Reagents and Kits for Capillary Ion Analysis

Description	Volume	Part No.
CIA-Pak® MBT Solution		WAT036991
ElectroOsmotic flow modifier for working electrolyte		
IonSelect Anion Package		WAT049385
Waters IonSelect high mobility anion electrolyte Waters IonSelect high mobility anion	1 L	
test solution	5 mL	
IonSelect Cation Package		WAT049386
Waters IonSelect low mobility cation electrolyte Waters IonSelect low mobility	1 L	
cation test solution	5 mL	
CIA-Pak UV Cat 1 Reagent		WAT054780
Alpha hydroxyisobutyric acid 15 mL bottle for CIA system	100 mL	
IonSelect OFM Hydroxide Concentrate		WAT049387
IonSelect 100 mM OFM hydroxide concentrate		
CIA-Pak HIBA		WAT054775
2-hydroxyisobutyric acid, 2 g		
MS Formic Acid Solution		700002341
Formic Acid	10 mL	
Triethylamine		100000110
Triethylamine used in Pico•Tag Chem Pkg Free AA analyses	100 mL	

Environmental Analysis

Description	Volume	Part No.
20 Pesticide Mix Standard		186006348
10 mg/mL of each pesticide: atrazin, atrazin-desethyl, cyanazin, sebuthylazin, simazin, terbuthylazin, hexazinone, chlortoluron, diuron, linuron, metobromuron, metoxuron, monolinuron, metolachlor, methamidophos, acephate, methomyl, dicrotophos, buprofenzin, fenpropimorph in acetonitrile	1 mL	
Carbamate Analysis Standards Set		186004626
Carbamate/Carbamoxyloxine Pesticides, Potable Water: This 2 mL flame-sealed ampule yields up to 2 Liters after dilution. Each standard contains all analytes below in the required range 15-150 µg/L. Aldicarb, Aldicarb sulfone, Aldicarb sulfoxide, Baygon, Carbaryl, Carbfuran, 3-Hydroxycarbonfuran, Methiocarb, Methomyl, Oxamyl (vydate) Carbamate Pesticides, Priority Pollutant: This 2 mL flame-sealed ampule yields up to 2 Liters after dilution. Standard contains a random selection from the list below in the required range 5-200 µg/L. Aldicarb, Aldicarb sulfone, Aldicarb sulfoxide, Baygon, Carbaryl, Carbfuran, Difuron, 3-Hydroxycarbonfuran, Methiocarb, Methomyl, Oxamyl (vydate), Propham	2 mL 2 mL	

The kits below should be used to assess the performance of the SFC or ACQUITY UPC 2m System. There are multi-component and individual component kits.

SFC and UPC² Standards

Description	Volume	Part No.
Waters Prep 15/30 SFC System Test Mix and Internal Standard		700005675
21 mg antipyrin, 45 mg ketoprofen, 15 mg sulfamethazine, 60 mL container 125 mg acetazolamide, 125 mL container	Solid Solid	
Waters Prep 100 SFC System Test Mix and Internal Standard		700005674
6250 mg antipyrin, 5625 mg ketoprofen, 625 mg fluorescein, 125 mL container 1250 mg sulfamethazine, 125 mL container	Solid Solid	
UPC ² Standard Mix		186006372
2 mg/mL each: 3-benzoylpyridine, cortisone, 4-nitroaniline, 4,4'-biphenol in methanol	1 mL	

Description	Volume	Part No.
UPC ² Gradient Standard		186006551
1 mg/mL coumarin, 1 mg/mL flavone, 2 mg/mL caffeine, 1 mg/mL thymine, 2 mg/mL prednisone in 2-propanol	1 mL	
UPC ² Caffeine Standard		186006614
1.0 mg/mL caffeine in 2-propanol	2 mL	
UPC ² Standard Ibuprofen		186006521
1.0 mg/mL ibuprofen in 2-propanol	2 mL	
UPC ² Standard Ketoprofen		186006522
1.0 mg/mL ketoprofen in 2-propanol	2 mL	
UPC ² Standard Flavone		186006523
1.0 mg/mL flavone in 2-propanol	2 mL	
UPC ² Standard Flurbiprofen		186006524
1.0 mg/mL flurbiprofen in 2-propanol	2 mL	

Ion Calibration Standards

Description	Part No.
125 mL Ammonia NH ₃	044
Ammonia as NH_3 in water, 1,000 mg/L	
125 mL Ammonia N	045
Ammonia as N in water, 1,000 mg/L	
125 mL Bromide	046
Bromide in water, 1,000 mg/L	
125 mL Chloride	047
Chloride in water, 1,000 mg/L	
125 mL Nitrate NO ₃	051-ERA
Nitrate as NO_3 in water, 1,000 mg/L	

Description	Part No.
125 mL Nitrate N	052-ERA
Nitrate as N in water, 1,000 mg/L	
125 mL Nitrite N	053
Nitrite as N in water, 1,000 mg/L	
125 mL Phosphate PO ₄	060
Phosphate as PO ₄ in water, 1,000 mg/L	
125 mL Phosphate P	061-ERA
Phosphate as P in water, 1,000 mg/L	
125 mL Sulphate	062
Sulphate in water, 1,000 mg/L	

Description	Part No.
125 mL Phosphorus	063
Phosphorous in water, 1,000 mg/L	
125 mL Silica	064
Silica in water, 1,000 mg/L	
125 mL Bromate	065
Bromate in water, 1,000 mg/L	
125 mL Chlorate	066
Chlorate in water, 1,000 mg/L	
125 mL Chlorite	067
Chlorite in water, 1,000 mg/L	
125 mL Perchlorate	068
Perchlorate in water, 1,000 mg/L	

QC Trace Metal Standards

Description	Part No.
125 mL Aluminum	011
Aluminum in nitric acid, 1,000 mg/L	
125 mL Antimony	012
Antimony in nitric acid, 1,000 mg/L	
125 mL Arsenic	013
Arsenic in nitric acid, 1,000 mg/L	
125 mL Barium	014
Barium in nitric acid, 1,000 mg/L	
125 mL Beryllium	015
Beryllium in nitric acid, 1,000 mg/L	
125 mL Boron	016
Boron in nitric acid, 1,000 mg/L	
125 mL Cadmium	017
Cadmium in nitric acid, 1,000 mg/L	
125 mL Calcium	018
Calcium in nitric acid, 1,000 mg/L	
125 mL Chromium VI	019
Chromium VI in water, 1,000 mg/L	
125 mL Chromium	020
Chromium in nitric acid, 1,000 mg/L	

Description	Part No.
125 mL Cobalt	021
Cobalt in nitric acid, 1,000 mg/L	
125 mL Copper	022
Copper in nitric acid, 1,000 mg/L	
125 mL Iron	023
Iron in nitric acid, 1,000 mg/L	
125 mL Lead	024
Lead in nitric acid, 1,000 mg/L	
125 mL Magnesium	025
Magnesium in nitric acid, 1,000 mg/L	
125 mL Manganese	026
Manganese in nitric acid, 1,000 mg/L	
125 mL Mercury	027
Mercury in nitric acid, 1,000 mg/L	
125 mL Molybdenum	028
Molybdenum in nitric acid, 1,000 mg/L	
125 mL Nickel	029
Nickel in nitric acid, 1,000 mg/L	
125 mL Potassium	030
Potassium in nitric acid, 1,000 mg/L	

Description	Part No.
125 mL Selenium	031
Selenium in nitric acid, 1,000 mg/L	
125 mL Silicon	032
Silicon in nitric acid, 1,000 mg/L	
125 mL Silver	033
Silver in nitric acid, 1,000 mg/L	
125 mL Sodium	034
Sodium in nitric acid, 1,000 mg/L	
125 mL Strontium	035
Strontium in nitric acid, 1,000 mg/L	
125 mL Thallium	036
Thalium in nitric acid, 1,000 mg/L	
125 mL Tin	037
Tin in hydrochloric acid, 1,000 mg/L	
125 mL Titanium	038
Titanium in hydrochloric acid, 1,000 mg/L	
125 mL Vanadium	039
Vanadium in nitric acid, 1,000 mg/L	
125 mL Zinc	040
Zinc in nitric acid, 1,000 mg/L	

Cation Analysis

Description	Part No.
125 mL Ammonium NH ₄	78102
Ammonium as NH_4 in water, 100 mg/L	
125 mL Ammonium N	78104
Ammonium as N in water, 100 mg/L	
125 mL Lithium	78110
Lithium in water, 100 mg/L	

Description	Part No.
125 mL Calcium	K10
Calcium in water, 1,000 mg/L	
125 mL Magnesium	K11
Magnesium in water, 1,000 mg/L	
125 mL Potassium	K12
Potassium in water, 1,000 mg/L	

Description	Part No.
125 mL Sodium	K13
Sodium in water, 1,000 mg/L	
125 mL Nitrite NO ₂	K15
Nitrite as NO_2 in water, 1,000 mg/L	

Suitability Standards

Description	Volume	Part No.
Neutral Suitability Standard		186006360
0.02 mg/L Thiourea, 0.25 mg/mL Naphthalene, 0.4 mg/mL Acenaphthalene in 50/50 acetonitrile/water	1 mL	
Reversed-Phase Suitability Standard		186006363
0.016 mg/mL uracil, 0.02 mg/mL butyl- paraben, 0.06 mg/mL naphthalene, 0.4 mg/mL propanol, 0.34 mgm/L dipropyl phthalate, 0.2 mg/mL acenaphthene, 0.1 mg/mL amitriptylene in 65/35 methanol/20mM K ₂ HPO ₄ buffer pH 7	1 mL	
Preparative Chromatography Standard		186006703
5 mg/mL Diclofenac sodium salt, 5 mg/mL Diphenhydramine hydrochloride, 5 mg/mL Flavone in DMSO	1 mL	
Dye Suitability Standard Kit		716000765
AutoPurification System Standard: (3) 2500 μg/mL thionin, 3000 μg/mL thioflavin, 2500 μg/mL crystal violet in 3/1 water/methanol	10 mL	

These pH Buffers are directly traceable to NIST SRMs, mercury free, guaranteed stable for at least one year after your receipt, and are supplied with a full certificate of analysis.

pH Buffers

Description	Volume	Part No.
pH 4 Liter		129
pH 4 Buffer	1 L	
pH 7 Liter		133
pH 7 Buffer	1 L	
pH 10 Liter		137
pH 10 Buffer	1 L	
pH 4 Pint		127
pH 4 Buffer	1 pint	
pH 7 Pint		131
pH 7 Buffer	1 pint	
pH 10 Pint		135
pH 10 Buffer	1 pint	

Trace Metal, Ion and Cation Calibration Standards can be used for primary calibration or to prepare second source calibration standards.

Find out more at www.waters.com/standards



Strengthen Your COMPLIANCE

It is becoming commonplace that global demands require work be performed to international recognized standards. Additionally, there is an increase in oversight and accreditation to these international standards. Waters Analytical Standards and Reagents make compliance simple and unambiguous. It also means that we become your standards and reagents compliance expert - included in the price of our products.

Laboratory Accreditation

Keeping up with laboratory accreditation requirements and audit practices is challenging for all organizations. Our products are manufactured to make that portion of laboratory accreditation and audits as easy as possible. Simply put, our products are challenged by auditors here at our plant, so they won't be at yours.

Product Traceability

All Waters Analytical Standards and Reagents are manufactured, traceable, and reproducible to exacting specifications accompanied by the highest level of accreditation and documentation. Only documented traceable products should be used in your laboratories, and ours are traceable, precisely formulated, prepared and tested by analytical chemists for analytical chemists.

Manufacturing to a new level

Waters has opened a brand new manufacturing facility designed to meet the stringent requirements that are expected in the manufacture and release of products associated with this business. Our facility is custom built for compliance, and fully accredited to ISO and GLP standards.

Analytical Standards and Reagents Accreditation

ISO 9001:2008 Registration

ISO/IEC 17025:2005 General Requirements for the Competence of Testing and Calibration Laboratories

ISO Guide 34:2009 General
Requirements for the Competence
of Reference Material Producers

Waters
THE SCIENCE OF WHAT'S POSSIBLE.





Chemical Materials

Whether you need to determine molecular weights, identify new compounds, monitor the effects of process or manufacturing changes, or troubleshoot customer and product support applications, Waters Standards can help you enhance the efficiency and quality of your chemical analyses. These high-quality, proven standards provide confidence in the results of key applications of your UPLC, UPC, UPC, UPLC/MS/MS, UPLC/QTof MS, and preparative/purification systems.

- Select from a broad range of products and kits
- Enhance efficiency and quality of chemical analysis
- Enjoy the superior speed, accuracy, flexibility and robustness of Waters technology

Gel Permeation Chromatography (GPC) Standards

Part No.
WAT011590
WAT011592
WAT011594
WAT011596
WAT011598
WAT011600
WAT011602
WAT011604

Description	Part No.
Polystyrene Standard 110,000	WAT011606
Polystyrene, 110,000 Da, 5 g	
Polystyrene Standard 180,000	WAT011608
Polystyrene, 180,000 Da, 5 g	
Polystyrene Standard 430,000	WAT011612
Polystyrene, 430,000 Da, 5 g	
Polystyrene Standard 780,000	WAT011614
Polystyrene, 780,000 Da, 5 g	
Polystyrene Standard 1,300,000	WAT011616
Polystyrene, 1,300,000 Da, 1 g	
Polystyrene Standard 2,800,000	WAT011618
Polystyrene, 2,800,000 Da, 1 g	
Polystyrene Standard	WAT011620

3,600,000

Polystyrene, 3,600,000 Da, 1 g

Description	Part No.
Polystyrene Standard 4,300,000	WAT011622
Polystyrene, 4,300,000 Da, 1 g	
Polystyrene Standard 5,200,000	WAT011624
Polystyrene, 5,200,000 Da, 1 g	
Polystyrene Standard 6,200,000	WAT011626
Polystyrene, 6,200,000 Da, 1 g	
Polystyrene Standard 8,400,000	WAT011628
Polystyrene, 8,400,000 Da, 1 g	
Polystyrene Standard 20,000,000	WAT011630
Polystyrene, 20,000,000 Da, 1 g	

Gel Permeation Chromatography (GPC) Standards Kits

For non-aqueous molecular weight determination.

Description	Part No.
Polymethylmethacrylate Mid MW Standards Kit	WAT035706
0.5 g/vial polymethylmethacrylate at each molecular weight: 2,400, 9,500, 31,000, 52,000, 10,0000, 170,000, 270,000, 490,000, 730,000, 1000000	
Polymethylmethacrylate Low MW Standards Kit	WAT035707
0.5 g/vial polymethylmethacrylate at each molecular weight: 1,000, 1,700, 2,500, 3,500, 5,000, 7,000, 10,000, 13,000, 20,000, 30,000	
Polyisoprene Standards Kit	WAT035708
0.5 g/vial polyisoprene at each molecular weight: 1000, 3,000, 10,000, 30,000, 70,000, 100,000, 300,000, 500,000, 1,000,000, 3,000,000	
Polybutadiene Standards Kit	WAT035709
0.5 g/vial polybutadiene at each molecular weight: 1,000, 3,000, 7,000, 10,000, 30,000, 70,000, 100,000, 300,000, 700,000, 1,000,000	
Polystyrene ReadyCal Standards 4 mL Kit	WAT058930
A complete kit of ready-to-use polystyrene calibration standards. Kit contains thirty 4 mL autosampler vials which contain four polystyrene standards per vial. There are three separate molecular weight ranges in each kit, ten	

units of each range. Range is from 400 to 2,000,000 Da.

Description	Part No.
Polystyrene ReadyCal Standards 2 mL Kit	WAT058931
A complete kit of ready-to-use polystyrene calibration standards. Kit contains thirty 2 mL autosampler vials which contain four polystyrene standards per vial. There are three separate molecular weight ranges in each kit, ten units of each range. Range is from 400 to 2,000,000 Da.	
Polystyrene Kit Low-Mid MW	WAT011588
A complete kit of polystyrene at MW: 400, 530, and 950, 10 g each MW: 2,800, 6,400, 10,000, 17,000, 43,000, 110,000 and 180,000, 5 g each	
Polystyrene Kit Mid-High MW	WAT011610
A complete kit of polystyrene at MW: 430,000, 780,000, 5 g each MW: 1300,000, 2,800,000, 3,600,000, 4,300,000, 5,200,000, 6,200,000, 8,400,000, and 2,000,0000, 1 g each	
Polystyrene Kit SL-105	WAT034208
A complete kit of polystyrene at MW: 520, 950, 1,200, 1,800, 2,470, 3,770, 5,100, 7,600, 12,500, 17,000, 500 mg each	
Polystyrene Kit SM-105	WAT034209
A complete kit of polystyrene at MW: 1200, 3250, 10,200, 28,000, 68,000, 195,000, 490,000, 1,080,000, 1,750,000, and 2,750,000, 500 mg each	
Polystyrene Kit SH-75	WAT034210
A complete kit of polystyrene at MW: 450,000, 1,270,000, 2,300,000, 3,260,000, 4,340,000, 8,000,000, and 15,000,000, 500 mg each	

The kits below should be used to assess the performance of the SFC or ACQUITY UPC $^{2^{tot}}$ System. There are multi-component and individual component kits.

SFC and UPC² Standards

Description	Volume	Part No.
Waters Prep 15/30 SFC System Test Mix and Internal Standard		700005675
21 mg antipyrin, 45 mg ketoprofen, 15 mg sulfamethazine, 60 mL container 125 mg acetazolamide, 125 mL container	Solid Solid	
Waters Prep 100 SFC System Test Mix and Internal Standard		700005674
6250 mg antipyrin, 5625 mg ketoprofen, 625 mg fluorescein, 125 mL container 1250 mg sulfamethazine, 125 mL container	Solid Solid	
UPC ² Standard Mix		186006372
2 mg/mL each: 3-benzoylpyridine, cortisone, 4-nitroaniline, 4,4'-biphenol in methanol	1 mL	
UPC ² Gradient Standard		186006551
1 mg/mL coumarin, 1 mg/mL flavone, 2 mg/mL caffeine, 1 mg/mL thymine, 2 mg/mL prednisone in 2-propanol	1 mL	
UPC ² Caffeine Standard		186006614
1.0 mg/mL caffeine in 2-propanol	2 mL	
UPC ² Standard Ibuprofen		186006521
1.0 mg/mL ibuprofen in 2-propanol	2 mL	
UPC ² Standard Ketoprofen		186006522
1.0 mg/mL ketoprofen in 2-propanol	2 mL	
UPC ² Standard Flavone		186006523
1.0 mg/mL flavone in 2-propanol	2 mL	
UPC ² Standard Flurbiprofen		186006524
1.0 mg/mL flurbiprofen in 2-propanol	2 mL	



ReadyCal Kits contain a full molecular weight range of polystyrene standards. Additional kits and individual GPC or GFC standards are available.

Learn more at www.waters.com/standards

Suitability Standards

Description	Volume	Part No.
•	volume	186006360
Neutral Suitability Standard		180000300
0.02 mg/L Thiourea, 0.25 mg/mL		
Naphthalene, 0.4 mg/mL Acenaphthalene in 50/50 acetonitrile/water	1 mL	
Reversed-Phase Suitability Standard		186006363
•		100000303
0.016 mg/mL uracil, 0.02 mg/mL butyl- paraben, 0.06 mg/mL naphthalene, 0.4 mg/mL propanol, 0.34 mgm/L dipropyl phthalate, 0.2 mg/mL acenaphthene, 0.1 mg/mL amitriptylene in 65/35		
methanol/20mM K_2HPO_4 buffer pH 7	1 mL	
Preparative Chromatography Standard		186006703
5 mg/mL Diclofenac sodium salt, 5 mg/mL Diphenhydramine hydrochloride,		
5 mg/mL Flavone in DMSO	1 mL	
Dye Suitability Standard Kit		716000765
AutoPurification™ System Standard: (3) 2500 μg/mL thionin, 3000 μg/mL thioflavin, 2500 μg/mL crystal violet in		
3/1 water/methanol	10 mL	

These pH Buffers are directly traceable to NIST SRMs, mercury free, guaranteed stable for at least one year after your receipt, and are supplied with a full certificate of analysis.

pH Buffers

Description	Volume	Part No.
pH 4 Liter		129
pH 4 Buffer	1 L	
pH 7 Liter		133
pH 7 Buffer	1 L	
pH 10 Liter		137
pH 10 Buffer	1 L	
pH 4 Pint		127
pH 4 Buffer	1 pint	
pH 7 Pint		131
pH 7 Buffer	1 pint	
pH 10 Pint		135
pH 10 Buffer	1 pint	

011	22, 48	06423, 48	18600328519, 36	18600661414, 47, 54
012	22, 48	06523, 48	18600328619, 36	18600670312, 21, 38,
013	22, 48	06623, 48	18600328719, 36	43, 49, 54
014	22, 48	06723, 48	18600382019, 36	18600684120, 37
015	22, 48	06823, 48	18600383617, 27, 33	18600684220, 37
016	22, 48	12714, 40, 43, 49, 54	18600383717, 27, 33	1860068615, 7
017	22, 48	12914, 40, 43, 49, 54	18600383817, 27, 33	1860068629
018	22, 48	13114, 40, 43, 49, 54	18600383917, 27, 33	1860068639
019	22, 48	13314, 40, 43, 49, 54	18600385622, 43	1860068649
020	22, 48	13514, 40, 43, 49, 54	18600385829, 43	1860068659
021	22, 48	13714, 40, 43, 49, 54	18600386319, 36	1860068669
022	22, 48	7810223, 49	18600386419, 36	18600686920, 37
023	22, 48	7810423, 49	18600413520, 37	18600687020, 37
024	22, 48	7811023, 49	18600462621, 47	20500034317, 33
025	22, 48	10000011029, 38, 47	18600490020, 37	43000178317, 33
026	22, 48	17600123517, 27, 33	18600600622, 29, 43	7000008896
027	22, 48	17600127917, 33	18600600722, 29, 43	70000107612
028	22, 48	17600253422, 43	18600600822, 43	7000013327
029	22, 48	17600298317, 33	18600601022, 43	7000015935, 7
030	22, 48	186000384C17, 27, 33	18600601118, 36	7000015945
031	22, 48	18600068928, 34	18600601218, 36	7000016735, 7
032	22, 48	186001860 18, 20, 36, 37	18600634821, 47	70000171313
033	22, 48	18600186118, 36	18600634920, 37	70000189213
034	22, 48	18600212218, 36	18600635021, 38	7000021576
035	22, 48	18600212318, 36	18600635122, 43	70000228613
036	22, 48	1860023258, 18,	18600636012, 21,	70000230913
037	22, 48	19, 35, 36	38, 43, 49, 54	7000023395, 7
038	22, 48	18600232618, 35	18600636312, 21, 38,	7000023407
039	22, 48	18600232718, 35	43, 49, 54	7000023417, 29, 38, 47
040	22, 48	1860023288, 18, 35	18600637118, 36	70000238713
044	23, 48	18600232918, 35	18600637214, 47, 54	70000239011
045	23, 48	18600233018, 35	18600651820, 37	70000239111
046	23, 48	18600233320, 37	18600651920, 37	7000024568
047	23, 48	18600233718, 35	18600652114, 47, 54	70000264211
051-ERA	23, 48	18600233818, 35	18600652214, 47, 54	7000026466
052-ERA	23, 48	18600278020, 37	18600652314, 47, 54	70000266911
053	23, 48	18600281720, 37	18600652414, 47, 54	7000027286
060	23, 48	18600286519, 36	18600655021, 38	7000027416
062	23, 48	18600286619, 36	18600655114, 47, 54	7000027527
063	23, 48	18600319619, 36	18600655220, 37	70000275312

700002820	7	70000567414, 47, 54	WAT01160223, 53	WAT04768512
700002837	7	70000567514, 47, 54	WAT01160423, 53	WAT04938529, 38, 47
700002846	11	71600076512, 21,	WAT01160623, 53	WAT04938629, 38, 47
700002848	11	38, 43, 49, 54	WAT01160823, 53	WAT04938729, 38, 47
700002857	7	71600323017, 33	WAT01161024, 53	WAT05287517, 27, 34
700002877	10	74100029917, 33	WAT01161223, 53	WAT05288027, 34
700002878	10	K1023, 49	WAT01161423, 53	WAT05288517
700002881	10	K1123, 49	WAT01161623, 53	WAT05289027, 34
700003093	7	K1223, 49	WAT01161823, 53	WAT05289527, 34
700003105	5, 7	K1323, 49	WAT01162023, 53	WAT05439224, 40
700003233	12	K1523, 49	WAT01162223, 53	WAT05477529, 38, 47
700003276	8	WAT00736028, 35	WAT01162423, 53	WAT05478029, 38, 47
700003694	12	WAT00757117, 27, 33	WAT01162623, 53	WAT05893024, 53
700003741	11	WAT01087329, 39	WAT01162823, 53	WAT05893124, 53
700003742	11	WAT01094728, 35	WAT01163023, 53	WAT08418928, 39
700003870	11	WAT01095028, 35	WAT03420724, 40	WAT08419828, 39
700003917	8	WAT01096028, 34, 35	WAT03420824, 53	WAT08419928, 39
700004021	10	WAT01096528, 34, 35	WAT03420924, 53	WAT08428228, 39
700004027	9	WAT01098328, 34	WAT03421024, 53	WAT08428328, 39
700004045	6	WAT01098528, 34	WAT03570624, 53	WAT08510128, 29, 39
700004155	12	WAT01157224, 40	WAT03570724, 53	WAT08510328, 29, 39
700004156	12	WAT01157424, 40	WAT03570824, 53	WAT08511028, 29, 39
700004380	11	WAT01157624, 40	WAT03570924, 53	WAT08514028, 29, 39
700004595	5, 7	WAT01157824, 40	WAT03571124, 40	WAT08514228, 29, 39
700004729	8	WAT01158024, 40	WAT03571424, 40	WAT08514429, 39
700004749	8	WAT01158224, 40	WAT03688829, 43	WAT08546628, 39
700004768	8	WAT01158424, 40	WAT03688922, 43	WAT08810828, 34, 35
700005089	8	WAT01158624, 40	WAT03699129, 38, 47	WAT08811228, 35
700005094	13	WAT01158824, 53	WAT03956622, 43	WAT08811928, 34, 35
700005095	13	WAT01159023, 53	WAT04287612	WAT08812029, 39
700005331	5	WAT01159223, 53	WAT04287912	WAT088122 17, 27, 33, 34
700005471	5	WAT01159423, 53	WAT04288112	WAT08812328, 34, 35
700005472	5	WAT01159623, 53	WAT04288512	WAT08813128, 35
700005666	9	WAT01159823, 53	WAT04288712	WAT09168128, 35
700005669	9	WAT01160023, 53	WAT04763412	

A		Atmospheric Pressure Photoionization Test Soluti	on 7
AccQ•Fluor Reagent Kit	17 27 24	AutoPurification System Standard	12, 21, 38, 43, 49, 54
AccQ•Tag Chemistry Kit		AutoSpec Premier Standards Test Kit	10
AccQ•Tag Eluent A		D	
• •		В	
AccQ•Tag Eluent B		Barium	22, 48
AccQ•Tag Ultra Chemistry Kit		BEH125 SEC Protein Standard Mix	20, 37
AccQ•Tag Ultra Derivatization Kit		BEH200 SEC Protein Standard Mix	20, 37
AccQ•Tag Ultra Eluent A		BEH450 SEC Protein Standard Mix	20, 37
AccQ•Tag Ultra Eluent B		Beryllium	22, 48
AccQ•Tag Ultra Chemistry Kit		Beverage Analysis 5 Standards Solution	22, 43
ACQUITY UPLC AAA Application Kit	17, 33	Beverage Analysis Kit	22, 43
ACQUITY UPLC/UV and ELSD		Beverage Analysis Mobile-Phase Reagent	22, 43
Qualification Standards Kit		Beverage Analysis Standard Solid	22, 43
ACQUITY UPLC/UV Qualification Standards Kit	11	Beverage Analysis Wash Reagent	
ACQUITY UPLC/UV Qualification Standards Kit		Beverage High Caffeine Standard	
for 3 or 5 mm Flow Cell		Boron	
ACQUITY UPLC/UV Qualification		Bromate	
Standards Kit 25 mm Flow Cell		Bromide	
ACQUITY UPLC/UV Start Up Standards Kit			20, 10
ACQUITY MS Test Solution	7, 29, 38, 47	C	
ACQUITY/Quattro micro or Quattro Premier MS		Cadmium	22, 48
Start Up Solution Kit	6	Caffeine Standard	12,
Alliance with Fluorescence Qualification	12	Calcium	•
Standards Test Kit	12	Calibration/Tuning on Autospecs	, ,
Alliance with RI Detector	11	Carbamate Analysis Standards Set	
Qualification Standards Kit		Cation Analysis	
Alliance HT Qualification Standards Kit		Chlorate	
Aluminum	22, 48	Chloride	•
Amino Acid Standard	17 27 22 24	Chlorite	ŕ
(AccQ•Tag, Pico•Tag, AccQ•Tag Ultra)			23, 48
Ammonia N	,	Cholesterol Sample for Atmospheric Pressure Photoionization Sensitivity Test	7
Ammonia NH ₃		Chromium	
Ammonium N	,		
Ammonium NH ₄		Chromium VI	
Analgesic Mix Standard		CI Optimization on GCTs	
Antimony	22, 48	CI /FI Optimization and Tuning on GCTs	
AquaAnalysis System Start Up Solution		CIAP & MDT C. A	
Arsenic	22, 48	CIAPak MBT Solution	
Atmospheric Pressure Ionization Test Solution		CIAPak UV Cat 1 Reagent	
(ZMD, QTof I/II/API-US/Quattro micro, Platform II, Quattro II	,	Cobalt	•
Atmospheric Pressure Photoionization Sensitivity Standard	I 7	Conductivity Detector OQ/PQ Solution	13

Copper	22, 48	Ion Select Anion Package	20, 37
Cytochrome C Digestion Standard	18, 36	Ion Select Cation Package	20, 37
2		Ion Select OFM Hydroxide Concentrate	20, 37
D		IonSABRE Positive Ion Sensitivity Standard	7
Dextrans Standard	24, 40	Iron	22, 48
Dextran Calibration Ladder	20, 37		
Dye Suitability Standard Kit 12, 21,	38, 43, 49, 54	L	
E		Lead	•
		Lithium	
Electrochemical Detector PQ Solution		LCT Premier/XE Standards Test Kit	9
ELS Detector Standard	13	M	
F		Magnesium	22, 23, 49
FI Optimization on GCTs	9	Manganese	22, 48
Fluorescence Detector Standard	12	Mass Scale and Resolution Set Up Solution	8
		MassPREP ADH Digestion Standard	8, 18, 35
G		MassPREP Bovine Hemoglobin Standard	18, 35
GC Headspace Qualification Standards Kit	13	MassPREP BSA Digestion Standard	18, 35
GC Performance Standard	13	MassPREP Digestion Standard Kit	18, 35
GC Qualification Standards Kit	13	MassPREP Digestion Standard Mix 1	19, 36
GCT Premier (El, Cl , DCl , SIP) Standards Kit	10	MassPREP Digestion Standard Mix 2	19, 36
GCT Premier (FI, FD) Standards Kit	10	MassPREP E. Coli Digest Standard	19, 36
Gel Permeation Chromatography Standards	23, 53	MassPREP Enhancer (5 vials)	19, 36
Gel Filtration Chromatography Standards Kits	24, 40	MassPREP Enolase Digest with Phosphopeptides Mix	19, 36
Gel Permeation Chromatography Standards Kits	24, 40	MassPREP Enolase Digestion Standard	8, 18, 35
[Glu 1]-Fibrinopeptide B Standard	8	MassPREP Glycoanalysis Kit	20, 37
Glycan Performance Test Standard	20, 37	MassPREP OST Standard	20, 37
Gradient Valve Standard	12	MassPREP Peptide Mixture	18, 35
		MassPREP Phosphopeptide Enrichment Kit	19, 36
Н		MassPREP Phosphopeptide Sample Kit - Enolase	19, 36
HPLC HT with UV Standards Kit	11	MassPREP Phosphopeptide Standard Enolase	19, 36
HPLC UV Linearity Solutions for Analytical Flow Cells	12	MassPREP Phosphorylase b Standard	18, 35
HPLC with UV Standards Kit	11	MassPREP Protein Standard Mix	20, 37
Hi3 E.Coli Standard	18, 36	MassPREP Enolase Digestion Standard	18, 35
Hi3 Phos B Standard	18, 36	Mercury	22, 48
		MetID Small Molecule Standard Mix	21, 38
		Molybdenum	22, 48
IEX Anion Exchange Standard		MS Calibration Reference Standard - PEG 400	5, 7
IEX Cation Exchange Standard	20, 37	MS Calibration Solution, Atmospheric Pressure Ionization	5, 7
Intact mAb Mass Check Standard	20, 37	MS ETD Reagent Kit	9
Ion Calibration Standards	23, 48	MS ETD Standards Kit	9


MS Leucine Enkephalin Kit	8	Pico•Tag Chem Pkg Free AA	28, 39
MS Qualification Standards Kit	6	Pico•Tag Diluent	28, 29
MS Reference Standard/Horse Heart Myoglobin	8	Pico•Tag Eluent 1	28, 39
MS Reference Standard/Reverse Peptide Kit	8	Pico•Tag Eluent 2	28, 39
MS Set Up Solution, Atmospheric Pressure Ionization	ı5, 7	Pico•Tag Eluent A	28, 39
MS Solution - Daidzein (1 mg/L)	7	Pico•Tag Eluent B	28, 39
MS Solution - Daidzein (250 μg/L)	5, 7	Pico•Tag Reagent Kit	28, 39
N		PITC Reagent	28, 39
N		Polyacrylic Acid Standards Kit	24, 40
Neutral Suitability Standard		Polybutadiene Standards Kit	24, 53
Nickel		Polyethylene Glycol Standards Kit	24, 40
Nitrate N		Polyethylene Oxide Kit	24, 40
Nitrate NO ₃		Polyethylene Oxide Standard 160,000	24, 40
Nitrite N		Polyethylene Oxide Standard 24,000	24, 40
Nitrite NO ₂	23, 48	Polyethylene Oxide Standard 340,000	24, 40
Р		Polyethylene Oxide Standard 40,000	24, 40
Perchlorate	23.49	Polyethylene Oxide Standard 570,000	24, 40
Performance and Operation Solutions for MS		Polyethylene Oxide Standard 79,000	24, 40
Pesticide Mix Standard		Polyethylene Oxide Standard 850,000	24, 40
pH 10 Liter		Polyisoprene Standards Kit	24, 53
pH 10 Pint		Polymethylmethacrylate Low MW Standards Kit	24, 53
pH 4 Liter		Polymethylmethacrylate Mid MW Standards Kit	24, 53
pH 4 Pint		Polystyrene ReadyCal Standards 2 mL Kit	24, 53
pH 7 Liter		Polystyrene ReadyCal Standards 4 mL Kit	24, 53
pH 7 Pint		Polystyrene Standard 1,300,000	23, 53
Phosphate P		Polystyrene Standard 10,100	23, 53
Phosphate PO4		Polystyrene Standard 110,000	23,53
Phosphorus	23, 48	Polystyrene Standard 17,000	23, 53
PIC A Low UV Reagent		Polystyrene Standard 180,000	23, 53
PIC B5 Low UV Reagent	28, 39	Polystyrene Standard 2,800	23, 53
PIC B6 Low UV Reagent	28, 39	Polystyrene Standard 2,800,000	23, 53
PIC B7 Low UV Reagent	28, 39	Polystyrene Standard 20,000,000	23, 53
PIC B8 Low UV Reagent	28, 39	Polystyrene Standard 3,600,000	23, 53
PIC Reagent A	28, 39	Polystyrene Standard 4,300,000	23,53
PIC Reagent B5	28, 39	Polystyrene Standard 400	23, 53
PIC Reagent B6	28, 39	Polystyrene Standard 43,000	23, 53
PIC Reagent B7	28, 39	Polystyrene Standard 430,000	23, 53
PIC Reagent B8	28, 39	Polystyrene Standard 5,200,000	23, 53
PIC Reagent D4	28, 39	Polystyrene Standard 530	23,53
PIC Reagent Sample Kit (PIC A, B5, B6, B7, B8)	28, 39	Polystyrene Standard 6,200,000	23, 53

Polystyrene Standard 6,400	23,53	Suitability Standards	12, 38, 43, 49, 54
Polystyrene Standard 780,000	23, 53	T	
Polystyrene Standard 8,400,000	23, 53		
Polystyrene Standard 950	23, 53	Thallium	,
Potassium	22, 23, 49	Tin	
Preparative Chromatography Standard	12, 21	Titanium	22, 48
Pullulan Kit	24, 40	TQD/Quattro Premier/XE/Quattro micro Kit	6
		Two Cuvette Kit	12
Q		U	
QC Trace Metal Standards	22, 48	UPC ² Caffeine Standard	14 47 54
QTof Qualification Standards Kit	8	UPC ² Flavone Standard	
QTof Standards Kit without Bovine	8		
Quadrapole: Tuning, Gradient, Sensitivity,		UPC ² Flurbiprofen Standard	
and Precision Standard	5, 7	UPC ² Gradient Standard	
Quattro micro GC Standards Kit	10	UPC ² Ibuprofen Standard	
R		UPC ² Ketoprofen Standard	
	10.00.00	UPC ² Standard Mix	
RapiGest	18, 20, 36, 37	UPC ² Standards Kit	
Refractive Index Performance and	10	UPLC AAA H-Class Applications Kit	
Operation Solution		UV Detector Linearity Standards Kit	
Reserpine Kit		UV Detector Linearity Standards Kit for 3 mm Flow Ce	
Resolution check on all GC/MS products		UV Detector PQ Standard	
Reversed-Phase Suitability Standard12,	21, 38, 43, 49, 54	UV Wavelength Accuracy (erbium perchlorate)	
S		UV Wavelength Accuracy (uracil)	12
SAM I Reagent	29 39	V	
Selenium		Vanadium	22 38
Silica		- Carlotterin	
Silicon	,	W	
Silver	,	Waters Prep 100 SFC System Test Mix	
Six Cuvette Set		and Internal Standard	14, 47, 54
Sodium		Waters Prep 15/30 SFC System Test Mix	
Soft Drink Analysis 4 Component Standard	, ,	and Internal Standard	14, 47, 54
·	22, 43	X	
Soft Drink Analysis 5 Component Standard and Mobile Phase Kit	22 //3		
Soft Drink Mobile Phase for NovaPak Column Analysis		Xevo TQ Standards Kit	
Soft Drink Mobile Phase for SunFire Column Analysis		Xevo TQ-S Set Up Standard	
SQD Set Up Solution		Xevo TQ-S Calibration and Resolution Standards	
·		Xevo TQ-S Qualification Standards Kit	5
SQD Standards Kit		Z	
Strontium	,		22.40
Sulphate	23, 48	Zinc	22, 48



Thank You for Your Interest in Waters Analytical Standards and Reagents. We Appreciate Your Business.

There are multiple ways in which you can place your order with Waters:

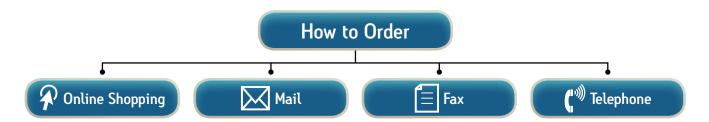

Online Shopping:

We encourage you to order online at www.waters.com/order

Waters offers a secure site, easy to use, available 24/7 for placing orders where you can:

- Check local stock availability
- Select multiple shipping options
- Upload directly excel sheet (CSV. format)
- Save time
- Save money by selecting promotional offers
- Get updated shipment information
- Track your business
- Save favorites and many more features...

By mail, fax, or telephone:

Contact your local Waters office, or visit www.waters.com/standards

SALES OFFICES:

Austria and European Export (Central South Eastern Europe, CIS and Middle East) 43 1 877 18 07

Belgium and Luxembourg 32 2 726 1000

Brazil 55 11 4134 3788

Canada 1 800 252 4752

China 86 21 6156 2666

Czech Republic 420 2 617 11384

Denmark 45 46 59 8080

Finland 358 9 5659 6288

France 33 1 30 48 72 00

Germany 49 6196 400 600

Hong Kong 852 2964 1800

Hungary 36 1 350 5086

India 91 80 2837 1900

Ireland 353 1 448 1500

Italy 39 02 265 0983

Japan 81 3 3471 7191

Korea 82 2 6300 4800

Mexico 52 55 52 00 1860

The Netherlands 31 76 508 7200

Norway 47 6 384 6050

Poland 48 22 101 5900

Puerto Rico 1 787 747 8445

Russia/CIS 7 495 727 4490/ 290 9737

Singapore 65 6593 7100

Spain 34 93 600 9300

Sweden 46 8 555 115 00

Switzerland 41 56 676 7000

US 1 800 252 4752

Milford, MA 01757 U.S.A. T: 508 478 2000 F: 508 872 1990



This seal is an emblem of quality, representing Waters' commitment to delivering standards and reagent products that meet our exacting requirements, and ensure the success of our customers.

www.waters.com/standards

Waters

THE SCIENCE OF WHAT'S POSSIBLE.™

